


RNAマーカー (100-1000 bases)
RNA marker (100-1000 bases)

Z230905A

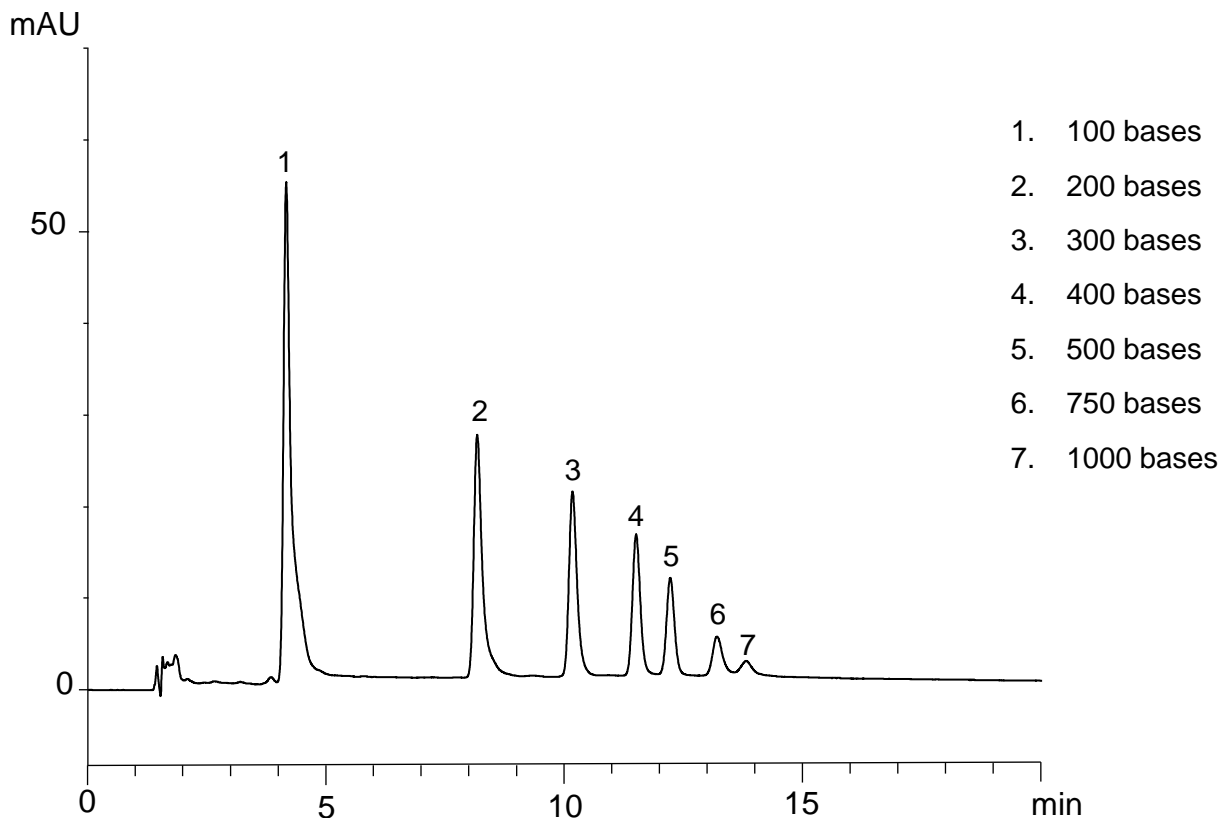

Column : Accura Triart Bio C18 (3 μ m, 30 nm)
100 X 2.1 mm I.D.

Eluent : A) 50 mM TEAA* (pH 7.0)
B) 50 mM TEAA* (pH 7.0)/acetonitrile (50/50)
18-27%B (0-20 min)

Flow rate : 0.2 mL/min

Temperature : 80°C

Detection : UV at 254 nm

Injection : 1 μ L (0.25 mg/mL)

Sample : Century™-Plus RNA Markers
(Thermo Fisher Scientific)

*triethylammonium acetate