
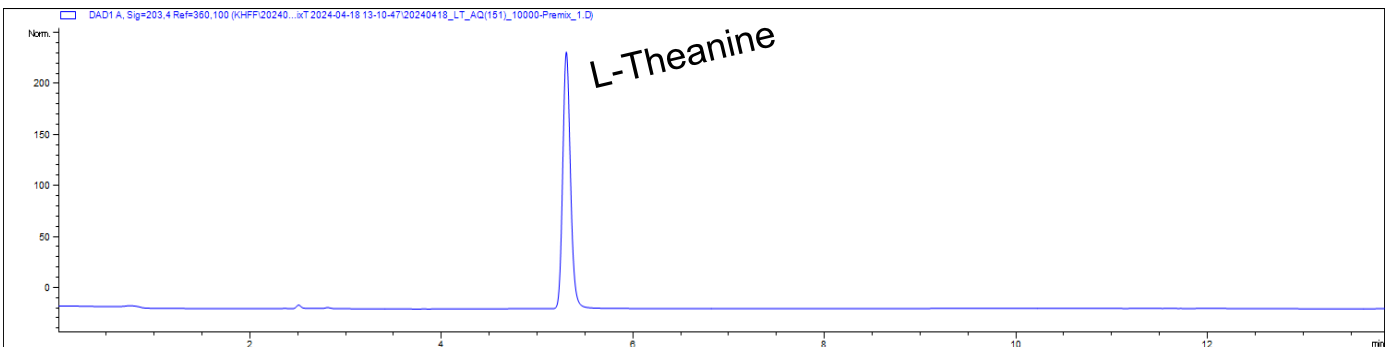


L-Theanine

테아닌

KR24050701

L-Theanine

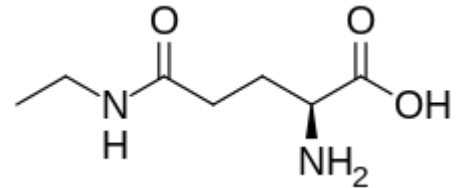

Column	: YMC-Pack ODS-AQ (5 $\mu$ m, 12nm) 250 x 4.6 mm I.D
Eluent	: A) 2% Acetonitrile in D.W B) Acetonitrile 0%B (0-10 min), 100%B (10.1-15.0min), 0%B (15.1-
Flow rate	: 20.0 min) 1.0 mL/min
Detection	: UV at 203 nm
Temperature	: 35 $^{\circ}$ C
Injection	: 10 $\mu$ L (L-Theanine : 100 $\mu$ L)