


脂溶性ビタミン
Oil-soluble vitamins

G910415B

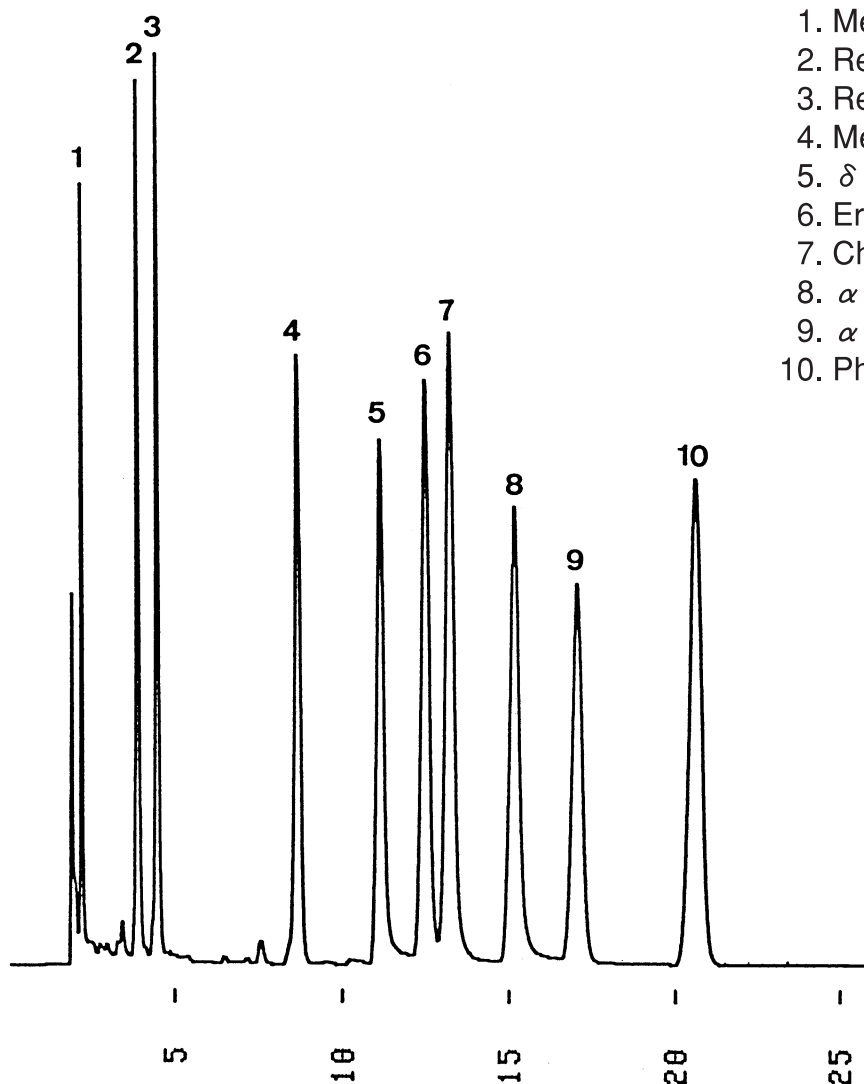

1. Menadione (Vitamin K₃)
2. Retinol (Vitamin A)
3. Retinol acetate
4. Menaquinone (Vitamin K₂)
5. δ -Tocopherol
6. Ergocalciferol (Vitamin D₂)
7. Cholecalciferol (Vitamin D₃)
8. α -Tocopherol (Vitamin E)
9. α -Tocopherol acetate
10. Phylloquinone (Vitamin K₁)

Column : YMC-Pack ODS-A (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 280 nm, 0.08 AUFS
Injection : 15 μ L (0.01~0.21 mg/mL)