
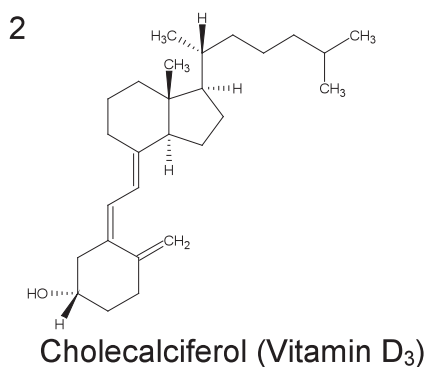
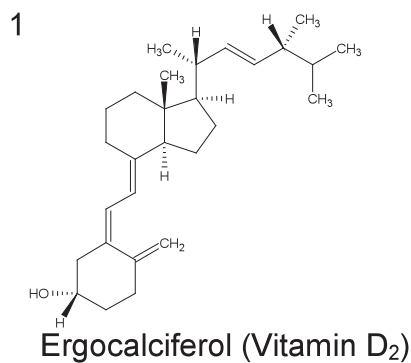


ビタミンD<sub>2</sub>とビタミンD<sub>3</sub>  
Vitamin D<sub>2</sub> and vitamin D<sub>3</sub>

P040819C

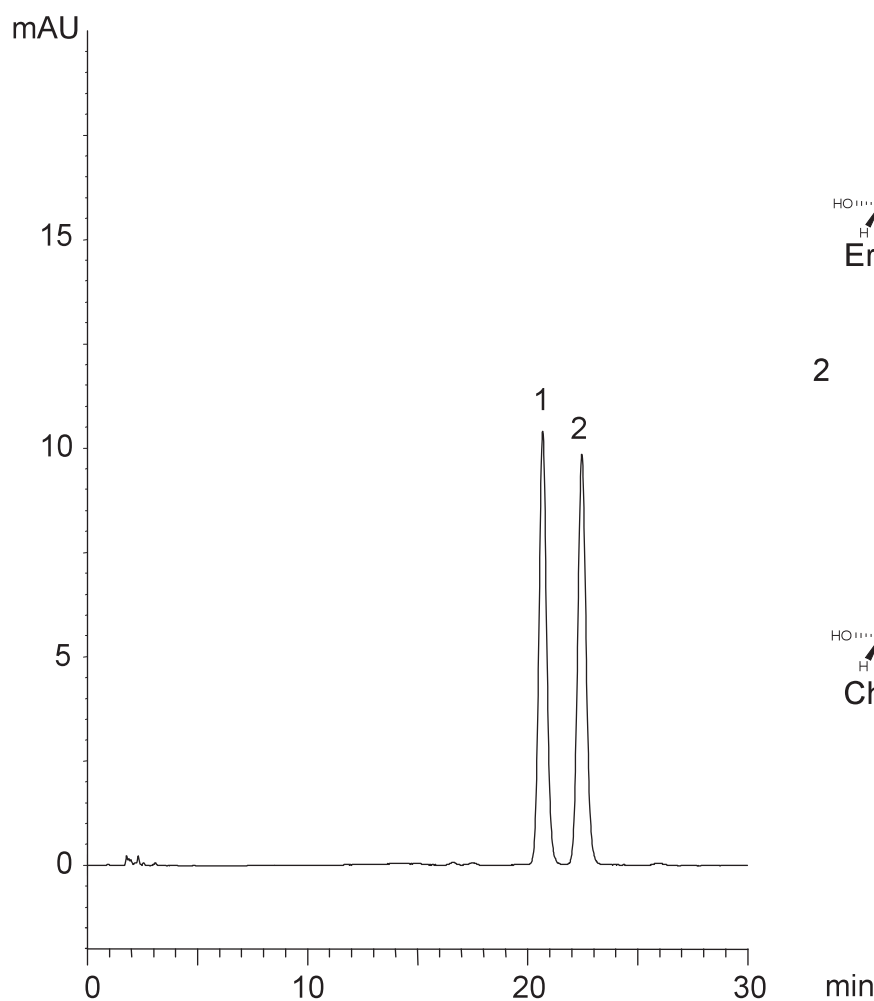

Column : YMC-Pack ODS-AL (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/water (95/5)  
Flow rate : 1.0 mL/min  
Temperature : 40  $^{\circ}$ C  
Detection : UV at 265 nm  
Injection : 10  $\mu$ L (0.01 mg/mL)