
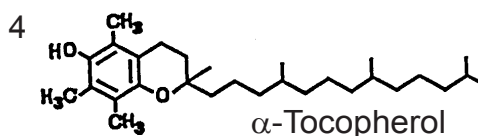
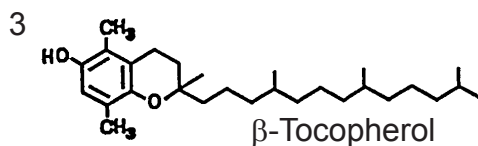
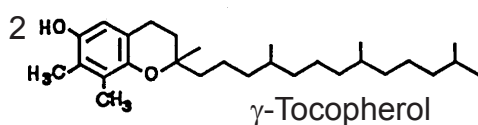
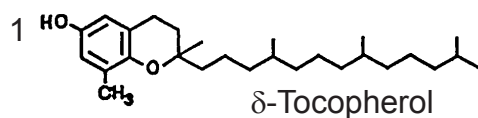


トコフェロール
Tocopherols

S921217A

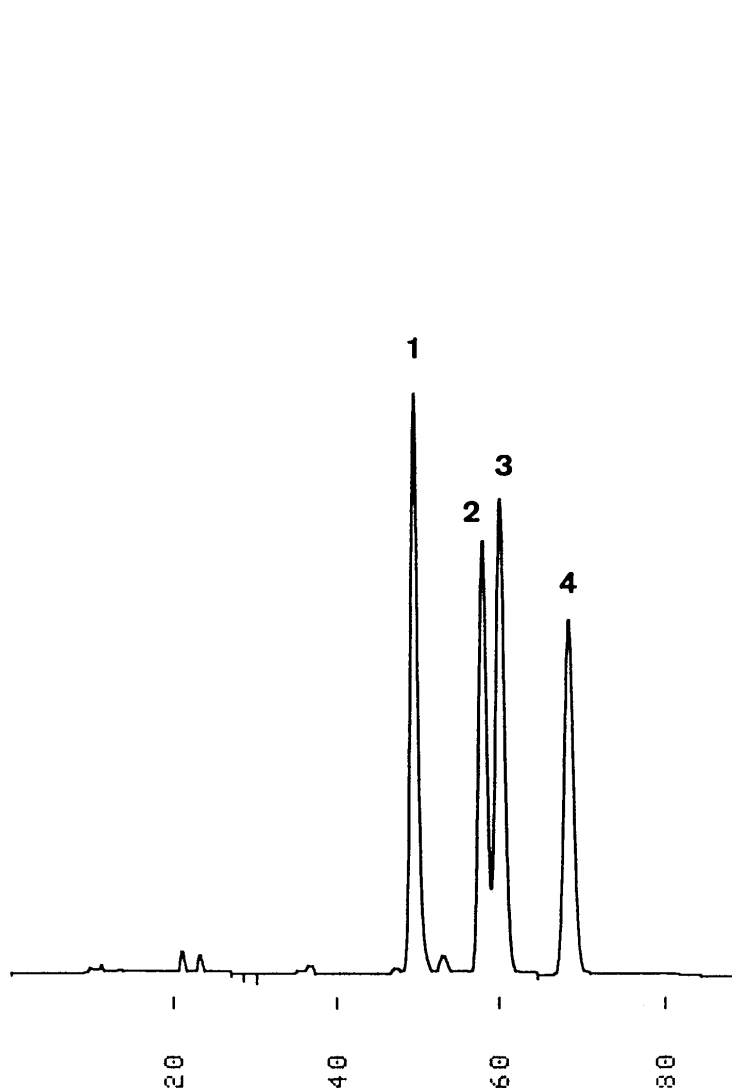

Column : YMC-Pack ODS-A (5 μm, 12nm)
250 × 4.6 mm I.D.

Eluent : 2-propanol/water (75/25)

Flow rate : 0.3 mL/min

Temperature : 30 °C

Detection : UV at 295 nm, 0.08 AUFS

Injection : 7 μL (0.25 mg/mL)

Ref.) Y.Satomura, M. Kimura and Y. Itokawa, J.Chromatogr., 625,
372-376 (1992)