
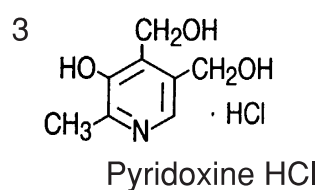
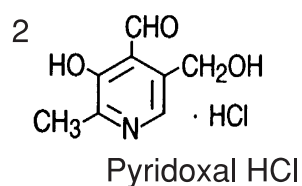
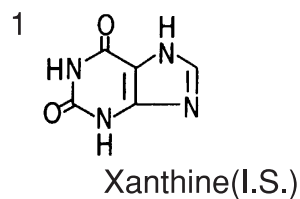


水溶性ビタミン  
Water-soluble vitamins

S960912B

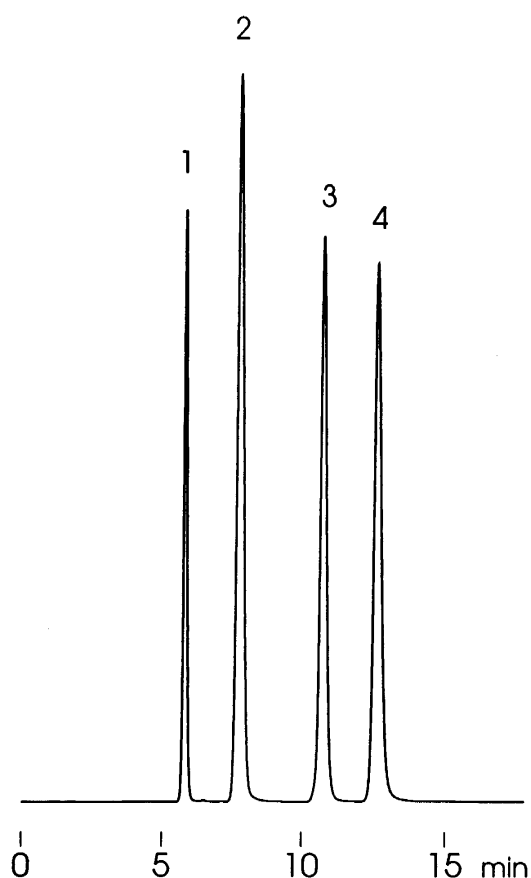

Column : YMC-Pack Pro C18 (5  $\mu$  m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : 20mM KH<sub>2</sub>PO<sub>4</sub> - K<sub>2</sub>HPO<sub>4</sub> (pH 6.9)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 254 nm, 0.13 AUFS  
Injection : 5  $\mu$ L (0.05~0.2 mg/mL)