


L-アスコルビン酸及びエリソルビン酸
L-Ascorbic acid and erythorbic acid

G910523N

1. L-Ascorbic acid (Vitamin C)
2. Erythorbic acid

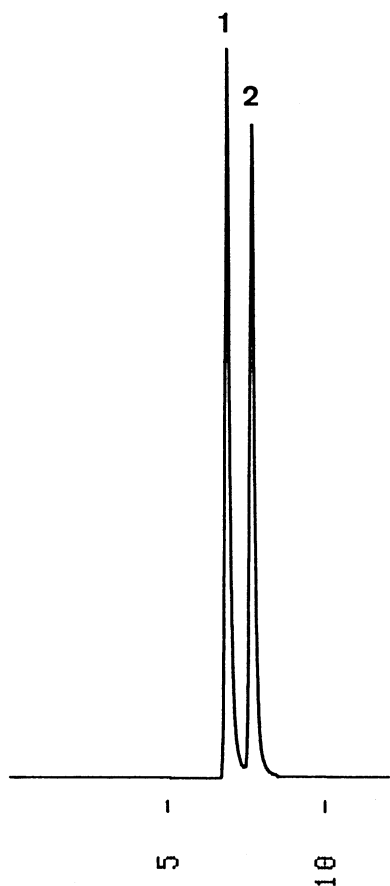

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 20mM CH₃COONH₄ in 5mM tetrabutylammonium bromide
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 250 nm, 0.08 AUFS
Injection : 5 μ L (0.05 mg/mL)