


オリゴ糖
Oligosaccharides

A970627F

1. Gal α 1 \rightarrow 6Glc α 1 \rightarrow 2 β Fruf (Raffinose)
2. [Gal α 1 \rightarrow 6]₂Glc α 1 \rightarrow 2 β Fruf (Stachyose)

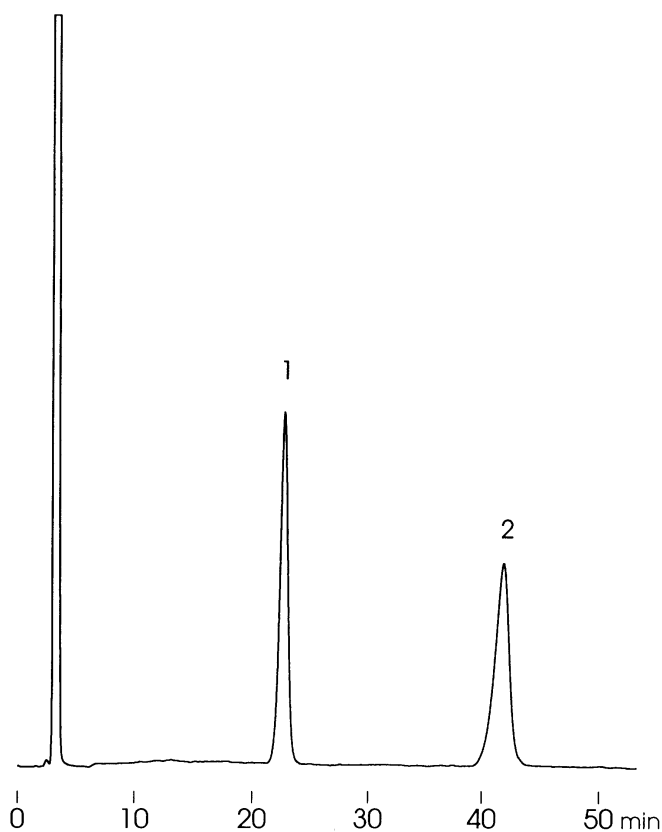

Column : YMC-Pack Polyamine II
250 × 4.6 mm I.D.
Eluent : acetonitrile/water (70/30)
Flow rate : 1.0 mL/min
Temperature : ambient (26 °C)
Detection : RI, 32 × 10⁻⁶ RIU/FS
Injection : 20 μL (4.0 mg/mL)