


オリゴ糖
Oligosaccharides

A970628H

1. Gal β 1 \rightarrow 4Glu α 1 \rightarrow 2 β Fruf
(Lactosyl fructoside)

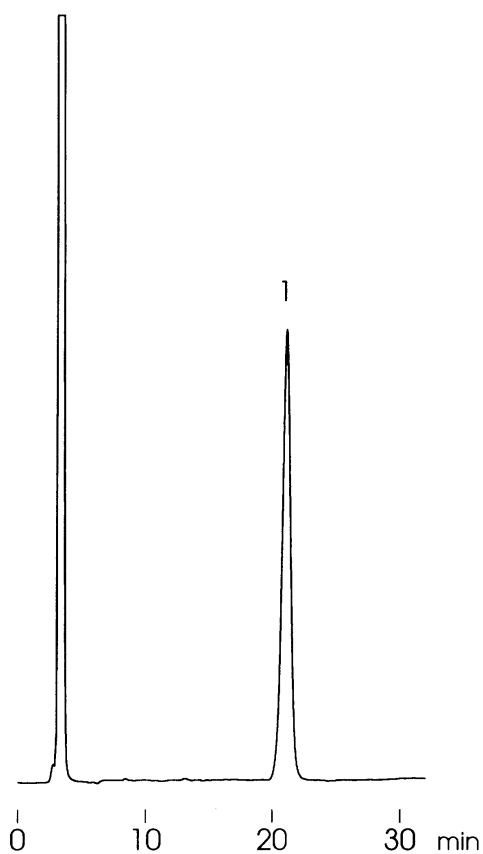

Column : YMC-Pack Polyamine II
250 × 4.6 mm I.D.
Eluent : acetonitrile/water (70/30)
Flow rate : 1.0 mL/min
Temperature : ambient (26 °C)
Detection : RI, 32 × 10⁻⁶ RIU/FS
Injection : 20 μL (4.0 mg/mL)