


マルトオリゴ糖(G1~G7)  
Maltooligosaccharides(G1~G7)

Y910501A

1. Glucose (G1)
2. Maltose (G2)
3. Maltotriose (G3)
4. Maltotetraose (G4)
5. Maltopentaose (G5)
6. Maltohexaose (G6)
7. Maltoheptaose (G7)

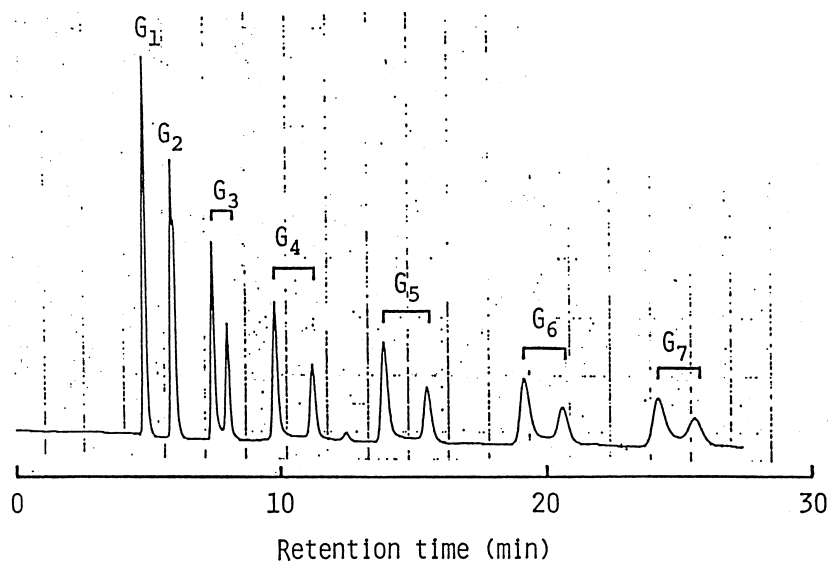

Courtesy of Dr.Kimura and Dr.Chiba, Hokkaido University, Japan

Column : YMC-Pack ODS-AQ (5  $\mu$  m, 12nm)  
250  $\times$  10.0 mm I.D.  
Eluent : water  
Flow rate : 3.0 mL/min  
Temperature : ambient (25  $^{\circ}$ C)  
Detection : RI