


マルトオリゴ糖(G11~G15)
Maltooligosaccharides(G11~G15)

M910604M

1. Maltoundecaose (G11)
2. Maltododecaose (G12)
3. Maltotridecaose (G13)
4. Maltotetradecaose (G14)
5. Maltopentadecaose (G15)

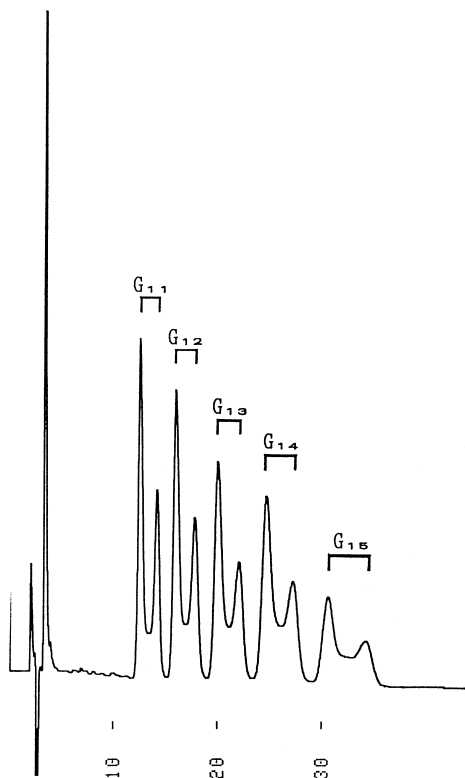

Column : YMC-Pack ODS-AQ (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : methanol/water (5/995)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : RI, 8×10^{-6} RIU/FS
Injection : 25 μ L (1 mg/mL)
Sample : Manufactured by Nakano Vinegar Co., Ltd.,
Nakano Central Research Institute