


糖ヌクレオチド
Sugar nucleotides

T911206D

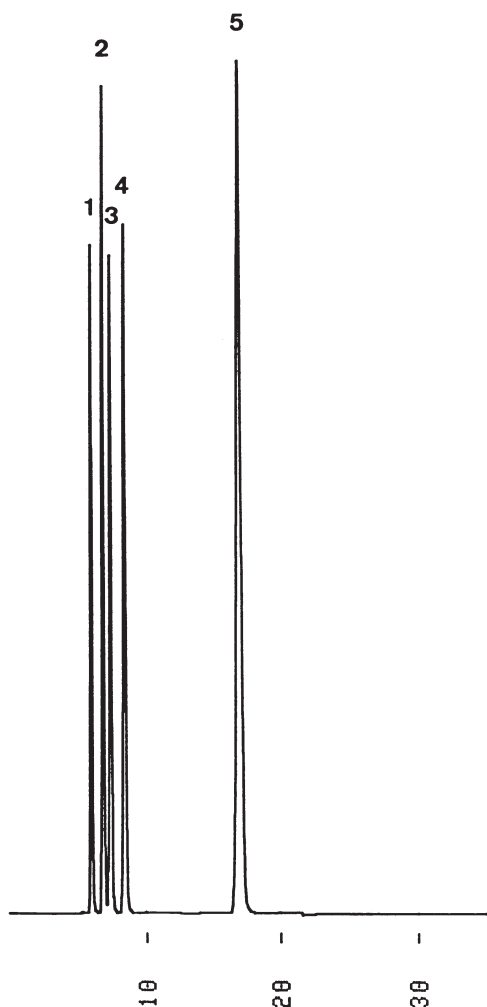

- 1.CDP-D-glucose
- 2.UDP-D-glucose
- 3.UDP-N-acetyl-D-glucosamine
- 4.GDP-D-mannose
- 5.ADP-D-glucose

Column : YMC-Pack Ph (5 μ m, 12nm)
250 \times 4.6 mm I.D.
Eluent : 100mM TEAA¹⁾ (pH 6.0)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 260 nm, 0.16 AUFS
Injection : 5 μ L (0.27~0.71 mg/mL)
1) TEAA: triethylamine-acetic acid