


アミノ酸
Amino acids

A000419B

1. L-Valine
2. L-Tyrosine
3. L-Phenylalanine
4. L-Tryptophan

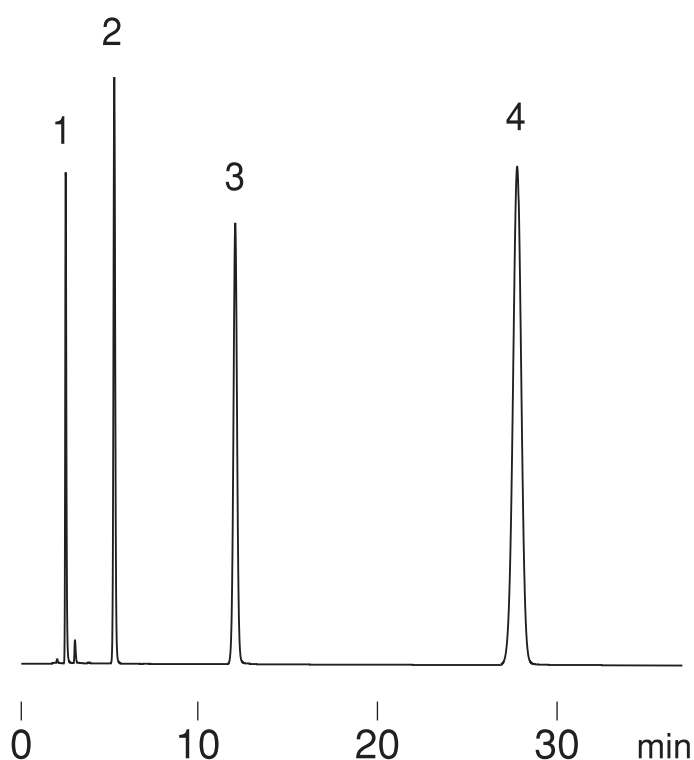

Column : Hydrosphere C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 20mM KH₂PO₄
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 210 nm, 0.16 AUFS
Injection : 5 μ L (0.05~1.6 mg/mL)