


PTH-アミノ酸
PTH-amino acids

K940331A

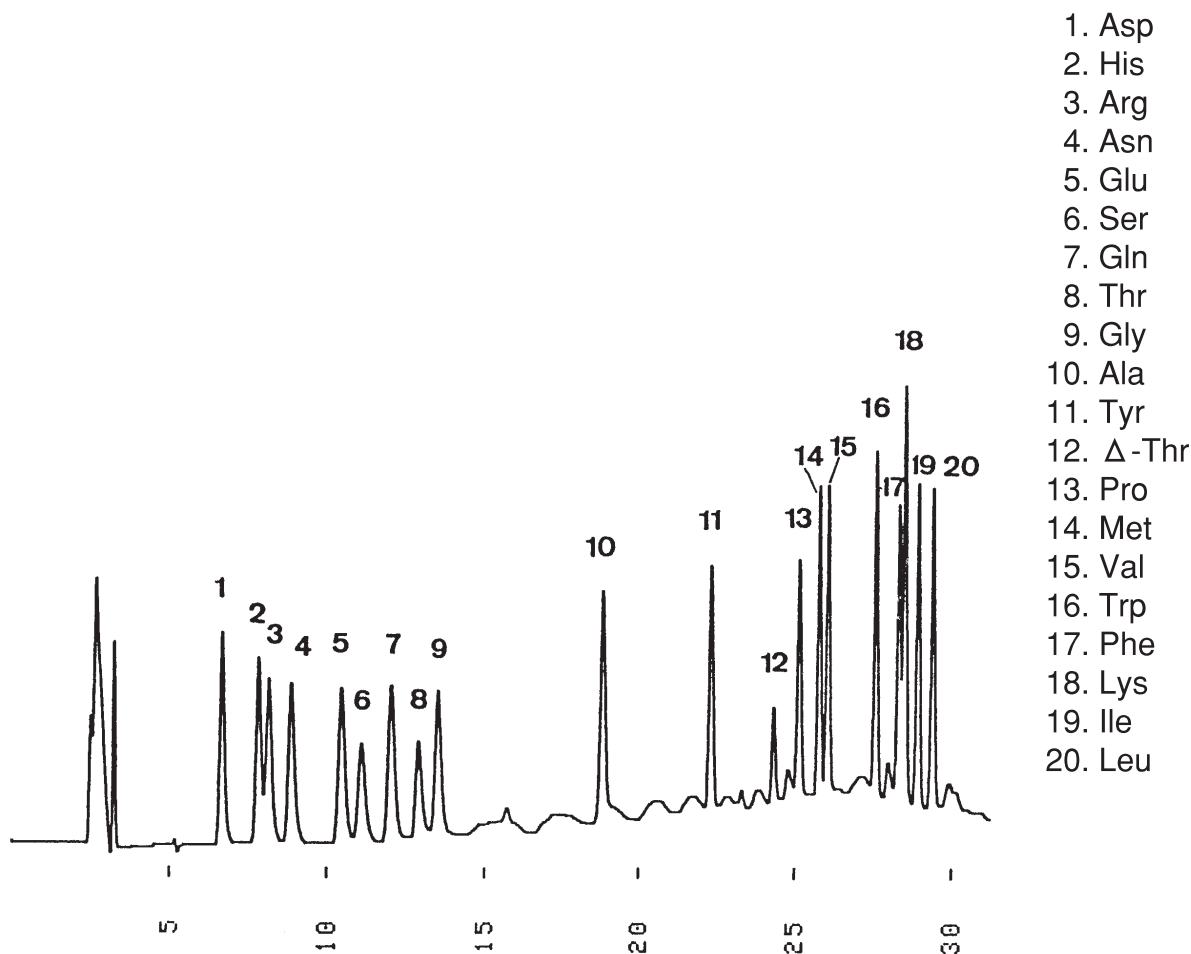

Column : J'sphere ODS-H80 (4 μ m, 8nm)
250 \times 2.0 mm I.D.

Eluent : A)acetonitrile/THF/300mM TEAA* (pH 5.1) (10/4/90)
B)acetonitrile/THF/300mM TEAA (pH 5.1) (75/4/25)
0-30%B (0-15min), 30-80%B (15-30min)

Flow rate : 0.2 mL/min

Temperature : 50 $^{\circ}$ C

Detection : UV at 270 nm, 0.08 AUFS

Injection : 30 μ L (25pmol/10 μ L)

*TEAA : triethylamine-acetic acid