


ペプチド
Peptides

A980827C

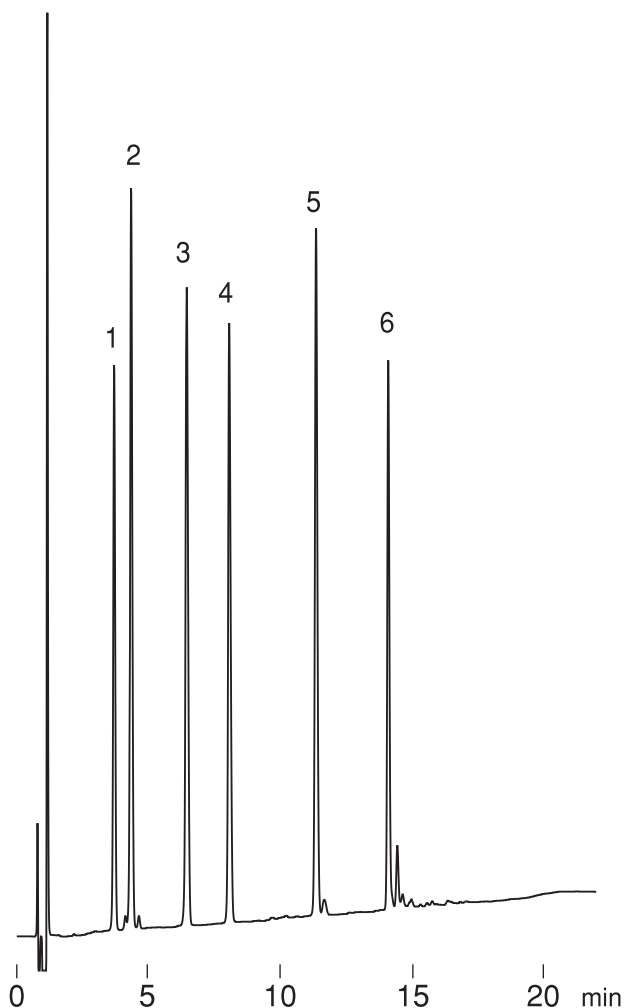

1. Oxytocin
2. Met-Enkephalin
3. Leu-Enkephalin
4. Angiotensin I
5. α -Mating factor
6. Insulin

Column : YMC-Pack Pro C18 (3 μ m, 12nm)
75 \times 4.6 mm I.D.

Eluent : A) water / TFA (100/0.1)
B) acetonitrile / TFA (100/0.1)
20-40%B (0-20min, linear)

Flow rate : 1.0 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 220 nm, 0.32 AUFS

Injection : 15 μ L (0.04~0.15 mg/mL)