


ペプチド
Peptides

A010821A

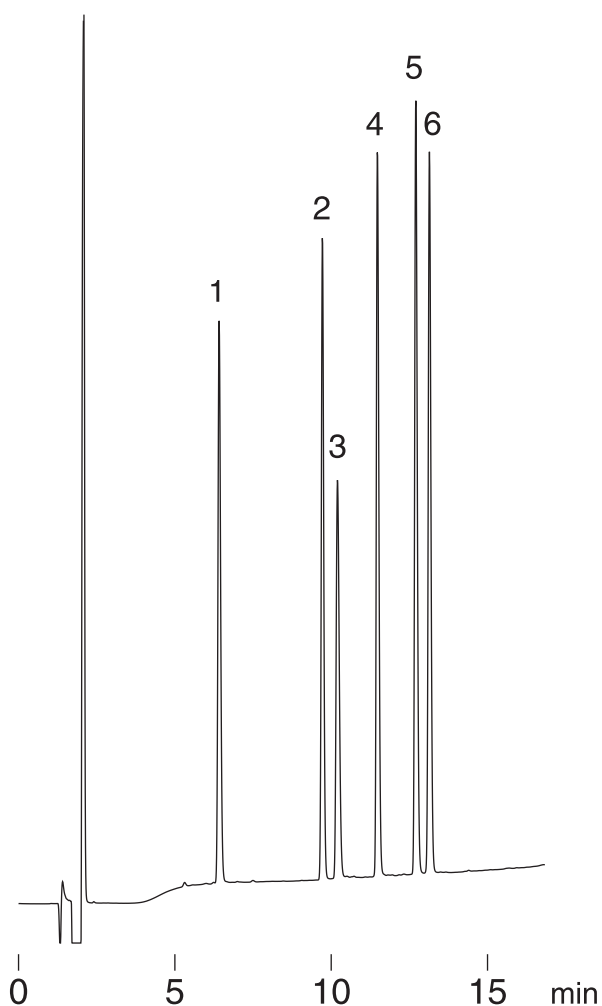

1. [Arg8]-Vasopressin
2. Oxytocin
3. Bradykinin
4. Angiotensin II
5. Leu-Enkephalin
6. Angiotensin I

Column : YMC-Pack Pro C18 RS (5 μ m, 8nm)
150 \times 4.6 mm I.D.

Eluent : A) water/TFA (100/0.1)
B) acetonitrile/TFA (100/0.1)
15-35%B(0-15 min)

Flow rate : 1.0 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 220 nm, 0.32 AUFS

Injection : 20 μ L (0.05 mg/mL)