


ペプチド  
Peptides

O940804C

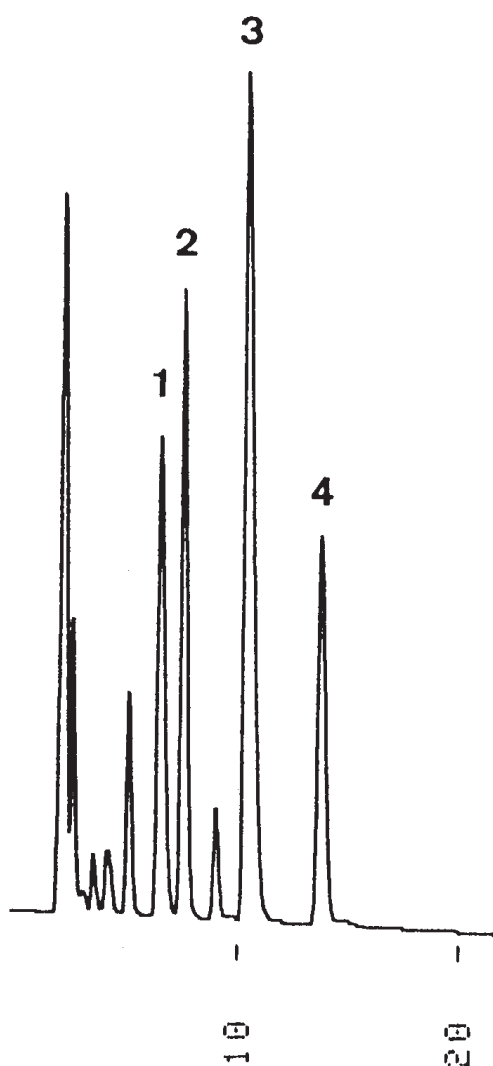

1. Bradykinin
2. Met-Enkephalin
3. Angiotensin III
4. Leu-Enkephalin

Column : J'sphere ODS-H80 (4  $\mu$ m, 8nm)  
150  $\times$  1.0 mm I.D.  
Eluent : acetonitrile/water/TFA (20/80/0.1)  
Flow rate : 0.05 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 220 nm, 0.08 AUFS  
Injection : 1  $\mu$ L (34~65  $\mu$ g/mL)