


ペプチド
Peptides

O940804D

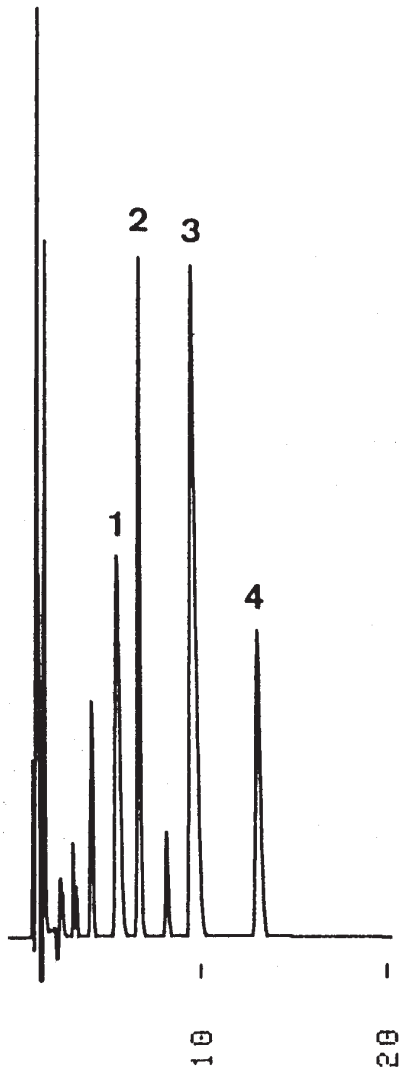

1. Bradykinin
2. Met-Enkephalin
3. Angiotensin III
4. Leu-Enkephalin

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : acetonitrile/water/TFA (20/80/0.1)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 220 nm, 0.32 AUFS
Injection : 20 μ L (34~65 μ g/mL)