


ペプチド (MW 556 - 3,465)
Peptides (MW 556 - 3,465)

R091009C

- | | |
|------------------------|------------|
| 1. Oxytocin | (MW 1,007) |
| 2. Met-Enkephalin | (MW 574) |
| 3. Leu-Enkephalin | (MW 556) |
| 4. Neurotensin | (MW 1,673) |
| 5. γ -Endorphin | (MW 1,859) |
| 6. β -Endorphin | (MW 3,465) |

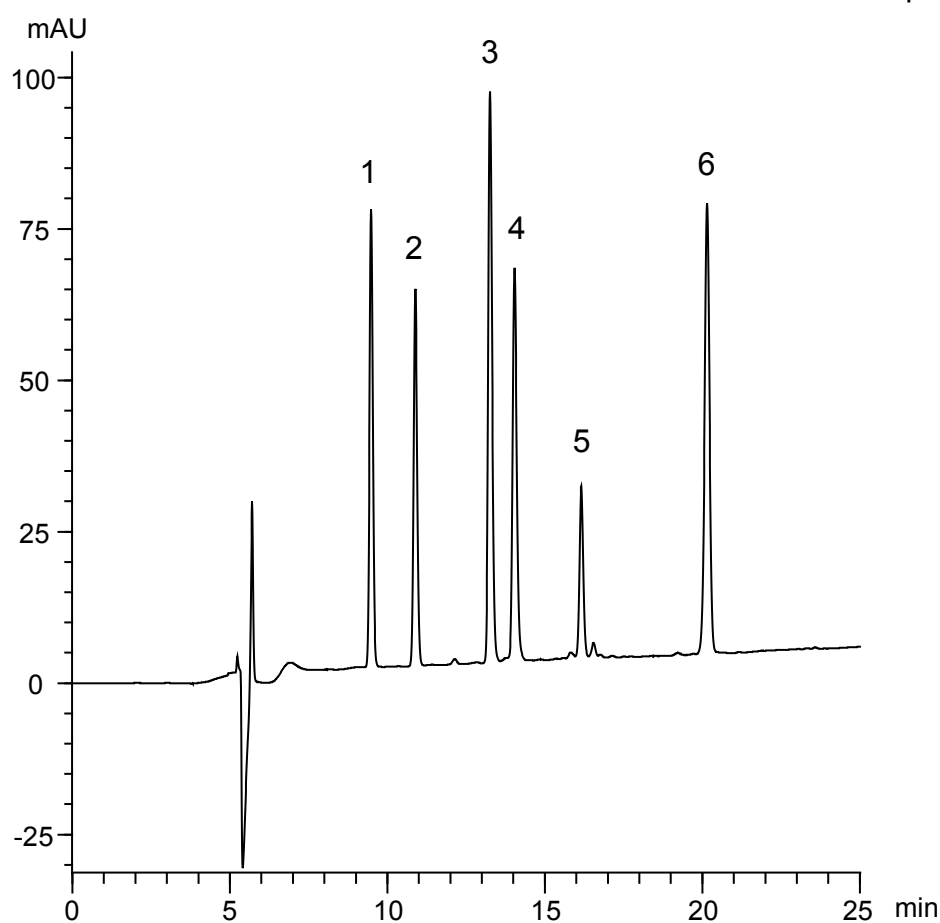

Column : YMC-Triart C18 (5 μ m, 12 nm)
150 X 2.0 mm I.D.

Eluent : A) water/TFA (100/0.1)
B) acetonitrile/TFA (100/0.1)
20-45%B (0-25 min)

Flow rate : 0.2 mL/min

Temperature : 37°C

Detection : UV at 220 nm

Injection : 2 μ L (0.075~0.25 mg/mL)