


タンパク質  
Proteins

N071016E

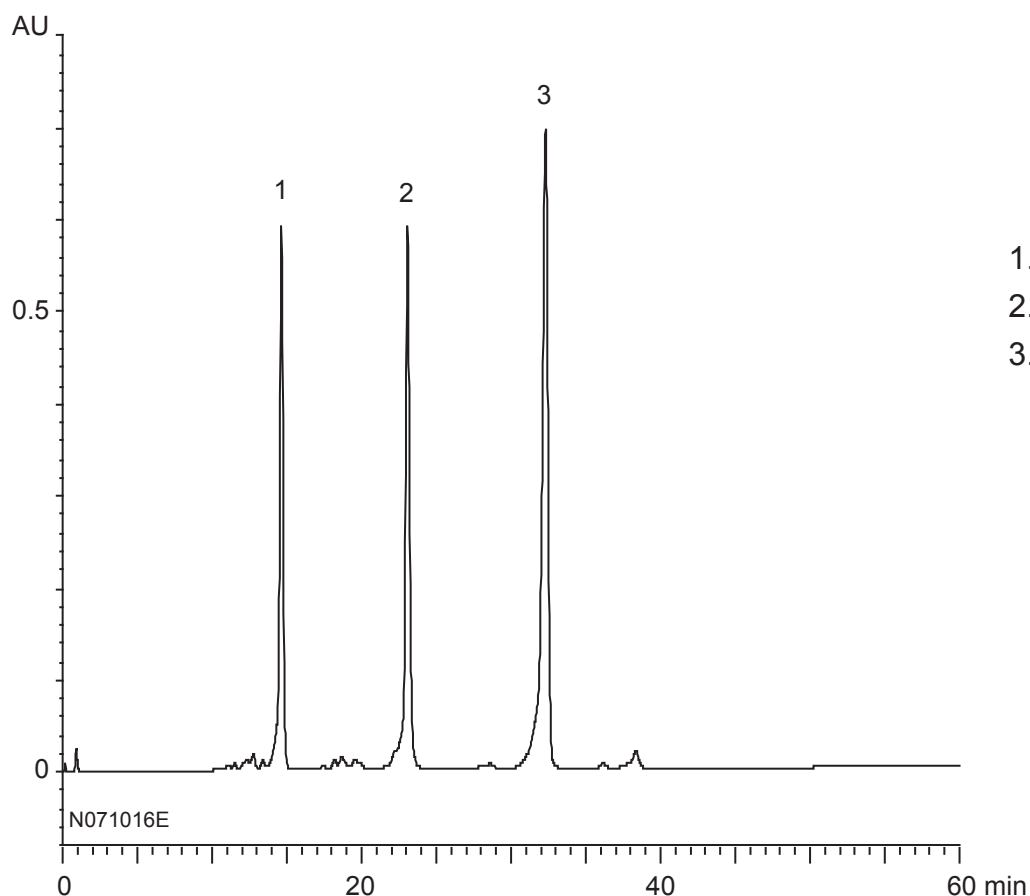

1. Ribonuclease A
2. Cytochrome c
3. Lysozyme

Column	: YMC-BioPro SP (5 $\mu$ m, 100 nm) 50 X 4.6 mm I.D.
Eluent	: A) 20 mM $\text{KH}_2\text{PO}_4$ - $\text{K}_2\text{HPO}_4$ (pH 6.8) B) 20 mM $\text{KH}_2\text{PO}_4$ - $\text{K}_2\text{HPO}_4$ (pH 6.8) containing 0.5 M NaCl 0-100%B (0-60 min)
Flow rate	: 0.5 mL/min
Temperature	: 25°C
Detection	: UV at 220 nm
Injection	: 20 $\mu$ L (0.5 mg/mL)