


タンパク質  
Proteins

N071023K

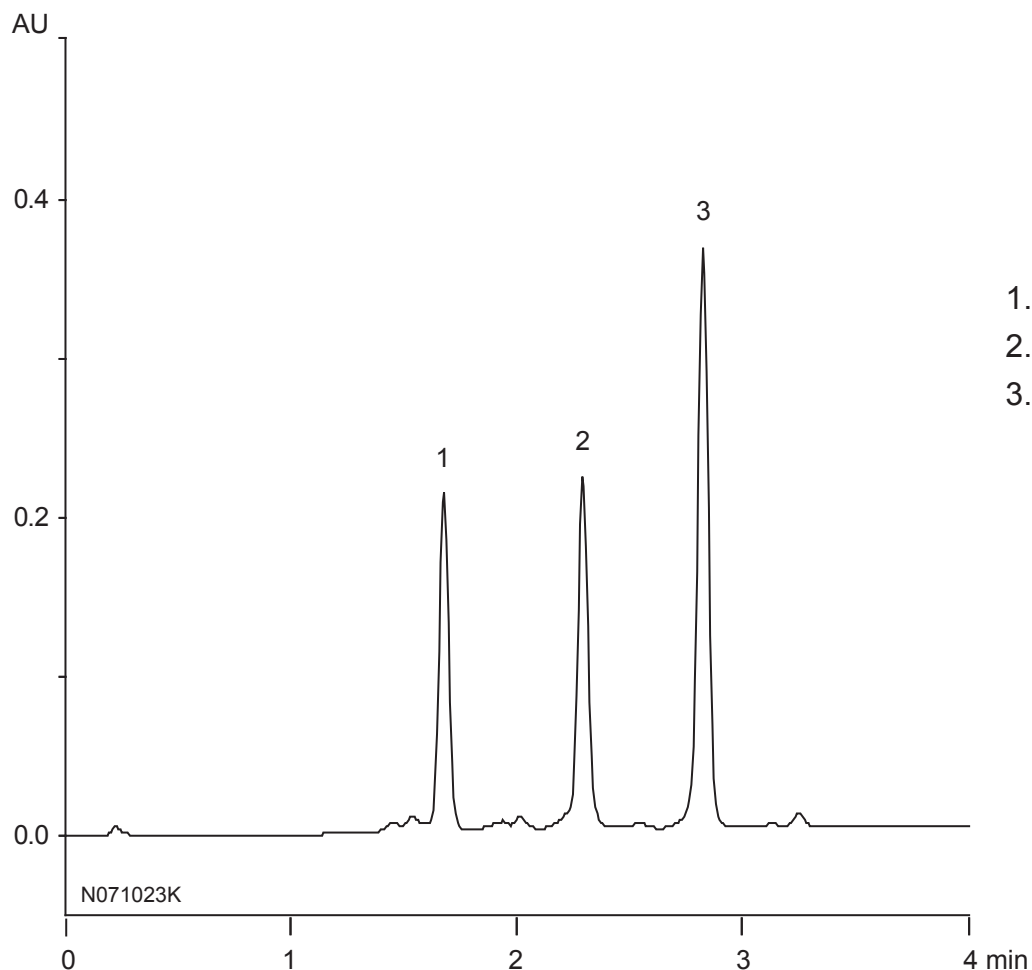

1. Ribonuclease A
2. Cytochrome c
3. Lysozyme

Column	: YMC-BioPro SP-F (5 μm) 30 X 4.6 mm I.D.
Eluent	: A) 20 mM KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 6.8) B) 20 mM KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 6.8) containing 0.5 M NaCl 0-100%B (0-4 min)
Flow rate	: 1.5 mL/min
Temperature	: 25°C
Detection	: UV at 220 nm
Injection	: 20 μL (0.1 mg/mL)