


ヒト血清中のタンパク質
Proteins in human serum

N080313E

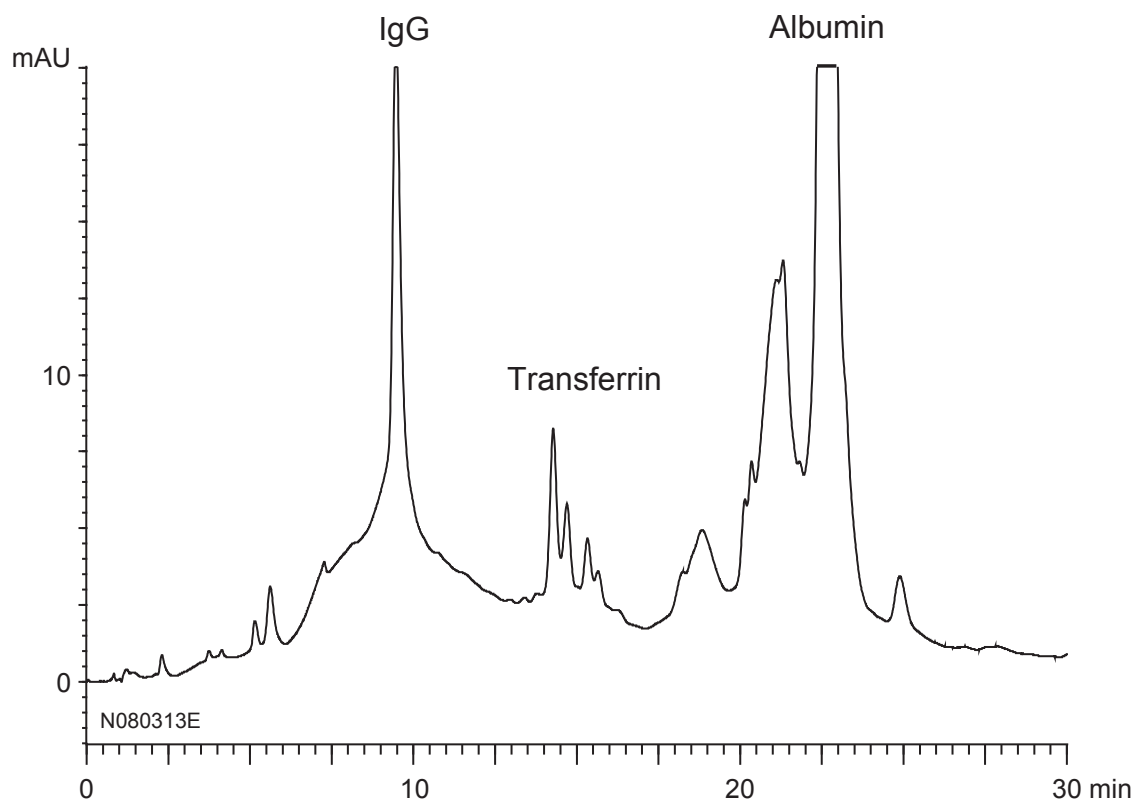

Column	: YMC-BioPro QA (5 μ m, 100 nm) 50 X 4.6 mm I.D.
Eluent	: A) 20 mM Tris-HCl (pH 8.6) B) 20 mM Tris-HCl (pH 8.6) containing 0.5 M NaCl 0-30%B (0-15 min), 30-100%B (15-30 min)
Flow rate	: 0.5 mL/min
Temperature	: 25°C
Detection	: UV at 280 nm
Injection	: 20 μ L
Sample	: Human serum (100 μ L/mL)