


モノクローナル抗体
Monoclonal antibody (MAb)

P080530A

MAb : Mouse monoclonal IgG1
(anti-human IgG4)

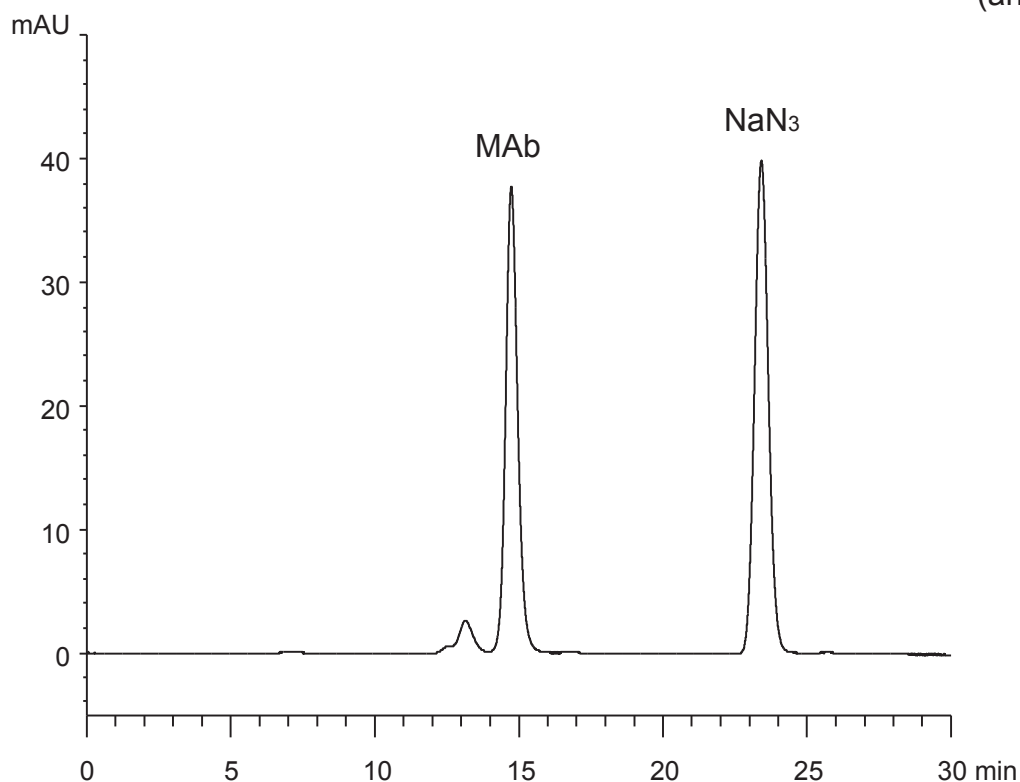

Column : YMC-Pack Diol-200
300 X 4.6 mm I.D.
Eluent : 0.1 M KH₂PO₄-K₂HPO₄ (pH 7.0)
Flow rate : 0.17 mL/min
Temperature : ambient (25°C)
Detection : UV at 220 nm
Injection : 10 µL
Sample : a commercially available mouse monoclonal IgG1 (0.05 mg/mL)
(purified by DEAE chromatography, containing NaN₃)