


モノクローナル抗体
Monoclonal antibody (MAb)

P080220A

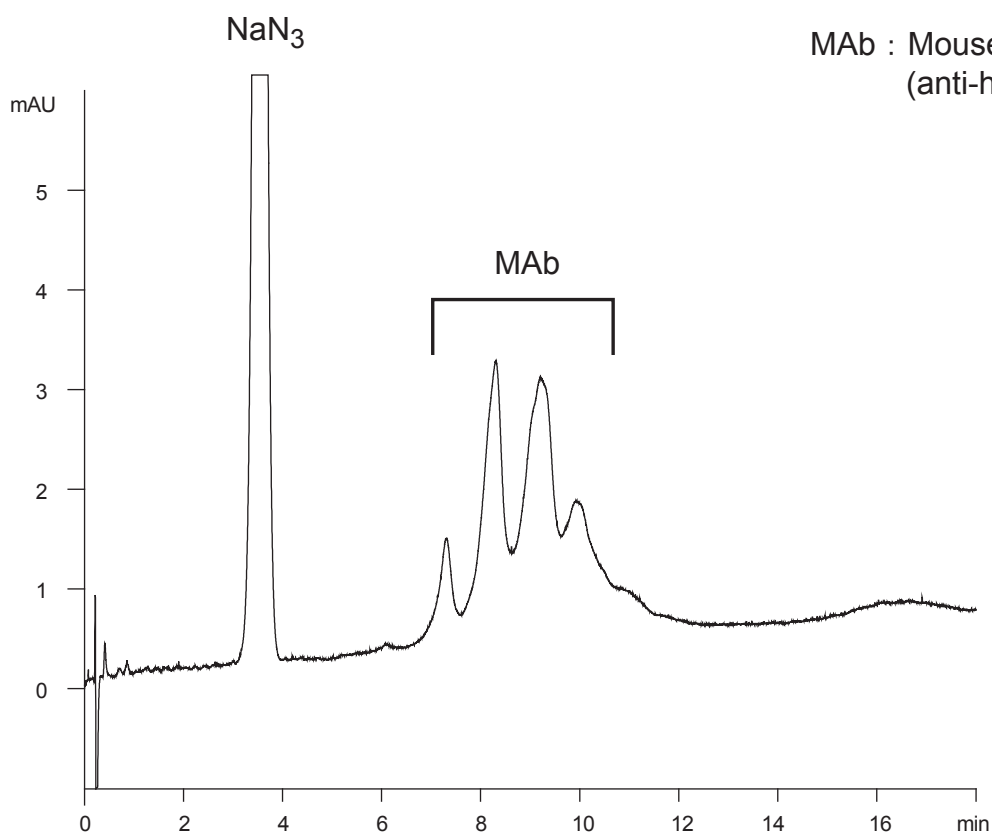

MAb : Mouse monoclonal IgG1
(anti-human IgG4)

Column : YMC-BioPro QA-F (5 μ m)
30 X 4.6 mm I.D.

Eluent : A) 20 mM Tris-HCl (pH 8.1)
B) 20 mM Tris-HCl (pH 8.1) containing 0.5 M NaCl
10-25%B (0-18 min)

Flow rate : 1.0 mL/min

Temperature : 25°C

Detection : UV at 220 nm

Injection : 10 μ L

Sample : a commercially available mouse monoclonal IgG1 (0.1 mg/mL)
(purified by DEAE chromatography, containing NaN₃)