


免疫グロブリン(IgG)および IgG フラグメント  
IgG, Fab and Fc fragments

R080708A

1. Rabbit IgG
2. Mouse IgG F(c) fragment
3. Mouse IgG Fab fragment

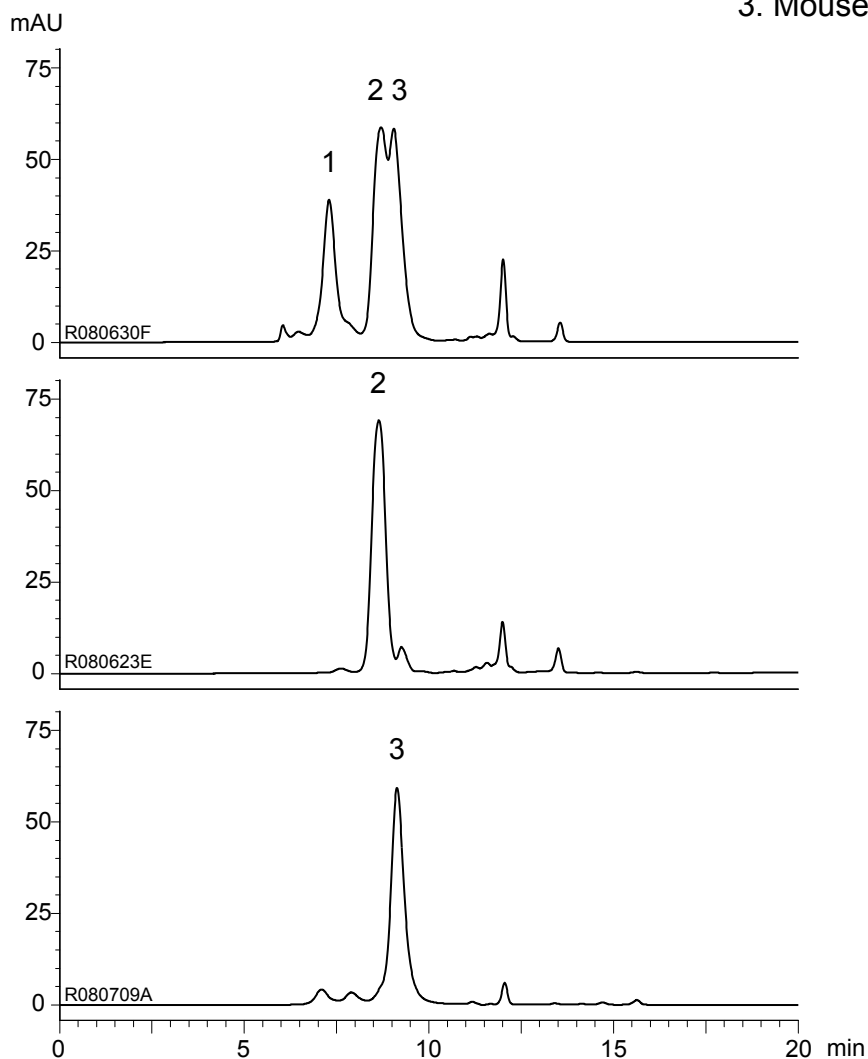

Column	: YMC-Pack Diol-200 300 X 8.0 mmI.D.
Eluent	: 0.1M KH <sub>2</sub> PO <sub>4</sub> -K <sub>2</sub> HPO <sub>4</sub> (pH 6.9) containing 0.2M NaCl
Flow rate	: 1.0 mL/min
Temperature	: ambient (27°C)
Detection	: UV at 220 nm
Injection	: 5 μL (0.4, 0.5 mg/mL)