


ヒト免疫グロブリン(IgM  $\kappa$  鎖, ミエローマタンパク質)  
Human IgM( $\kappa$ ), myeloma protein

G920108D

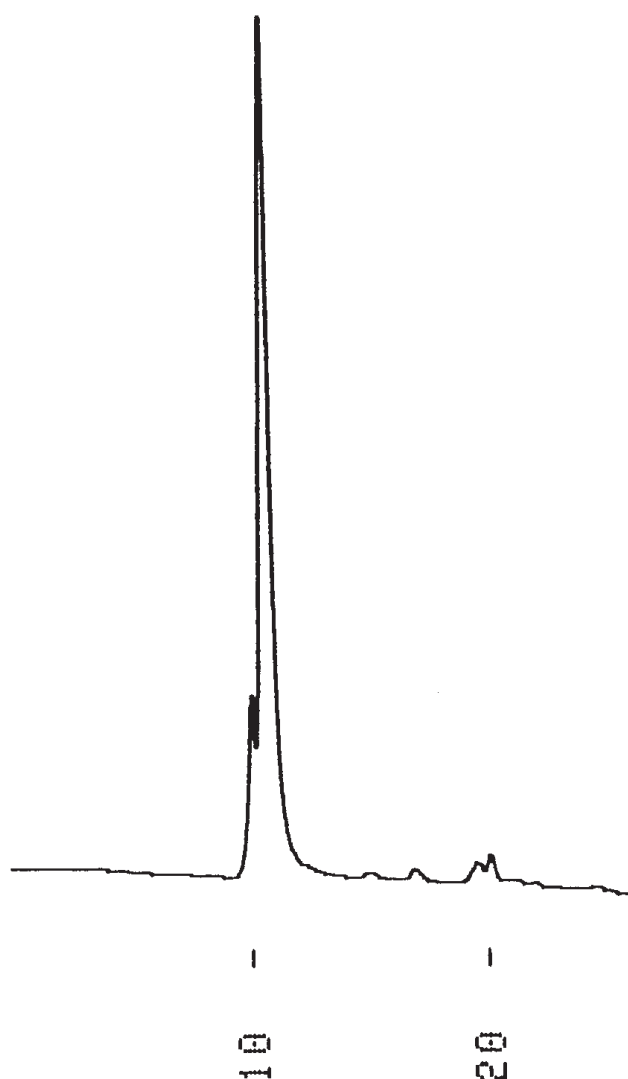

1. Normal human IgM ( $\kappa$ ), myeloma protein

Column	: YMC-Pack Diol-300 500 × 8.0 mm I.D.
Eluent	: 0.1M NaH <sub>2</sub> PO <sub>4</sub> - Na <sub>2</sub> HPO <sub>4</sub> (pH 6.8) containing 0.1M Na <sub>2</sub> SO <sub>4</sub>
Flow rate	: 1.0 mL/min
Temperature	: ambient (24 °C)
Detection	: UV at 280 nm, 0.04 AUFS
Injection	: 20 $\mu$ L (1.0 mg/mL)