


タンパク分子量マーカー

Proteins for molecular weight marker

G911214H

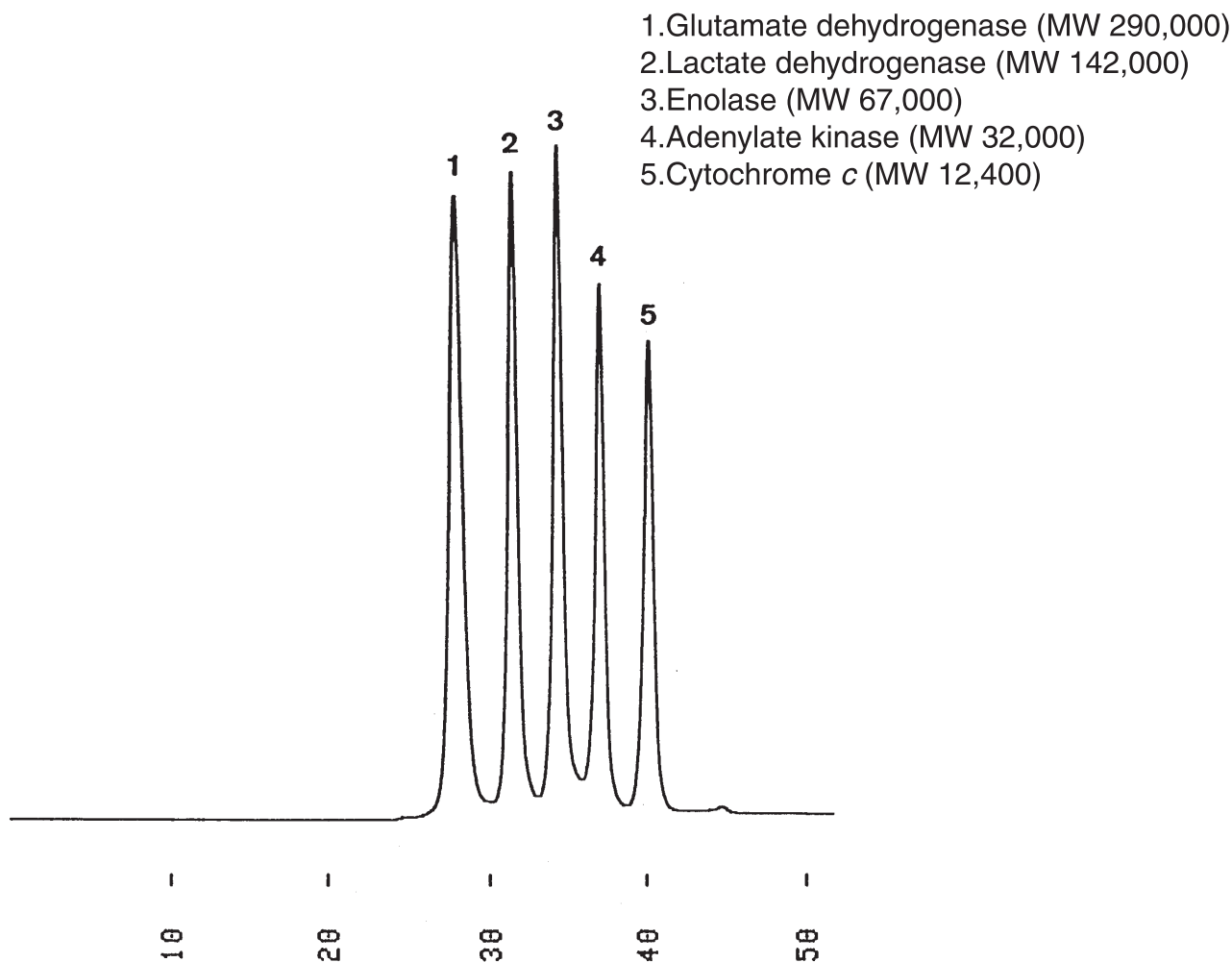

Column : YMC-Pack Diol-200 + Diol-200  
 500 × 8.0 mm I.D. + 300 × 8.0 mm I.D.  
 Eluent : 0.1M KH<sub>2</sub>PO<sub>4</sub>-K<sub>2</sub>HPO<sub>4</sub> (pH 7.0) containing 0.2M NaCl  
 Flow rate : 0.7 mL/min  
 Temperature : ambient (26 °C)  
 Detection : UV at 280 nm, 0.08 AUFS  
 Injection : 15 μL (100 μL/1 vial)  
 Sample : MW-Marker(HPLC)  
 manufactured by ORIENTAL YEAST CO.,LTD.