


不飽和脂肪酸メチルエステル
Unsaturated fatty acid methyl esters

K930323A

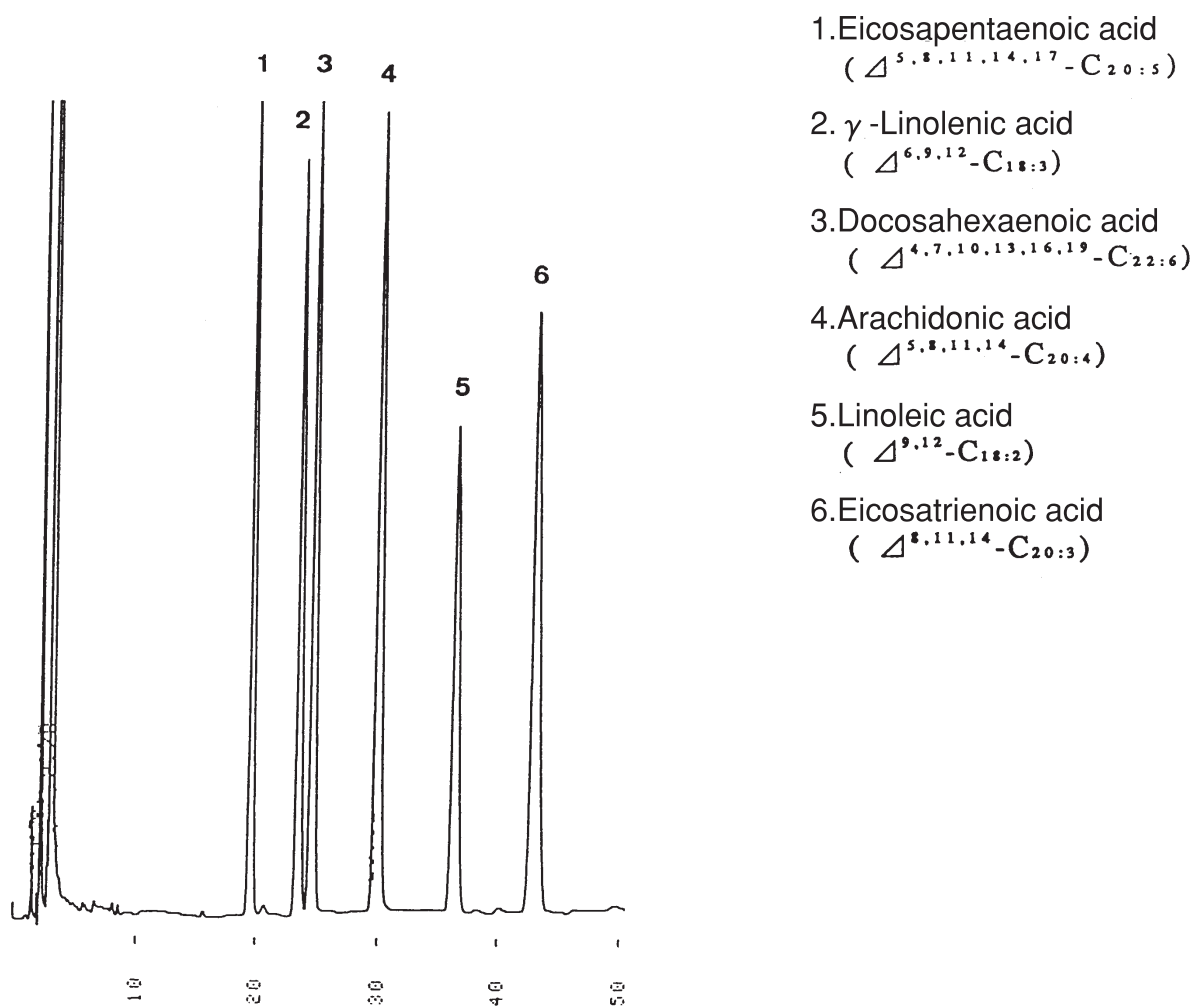

Courtesy of Dr. Yuji Saito, Tokyo Institute of Technology

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
250 \times 4.6 mm I.D.

Eluent : acetonitrile/water (81/19)

Flow rate : 1.5 mL/min

Temperature : ambient

Detection : UV at 190 nm, 0.128 AUFS

Injection : 20 μ L (0.165 mg/mL)