
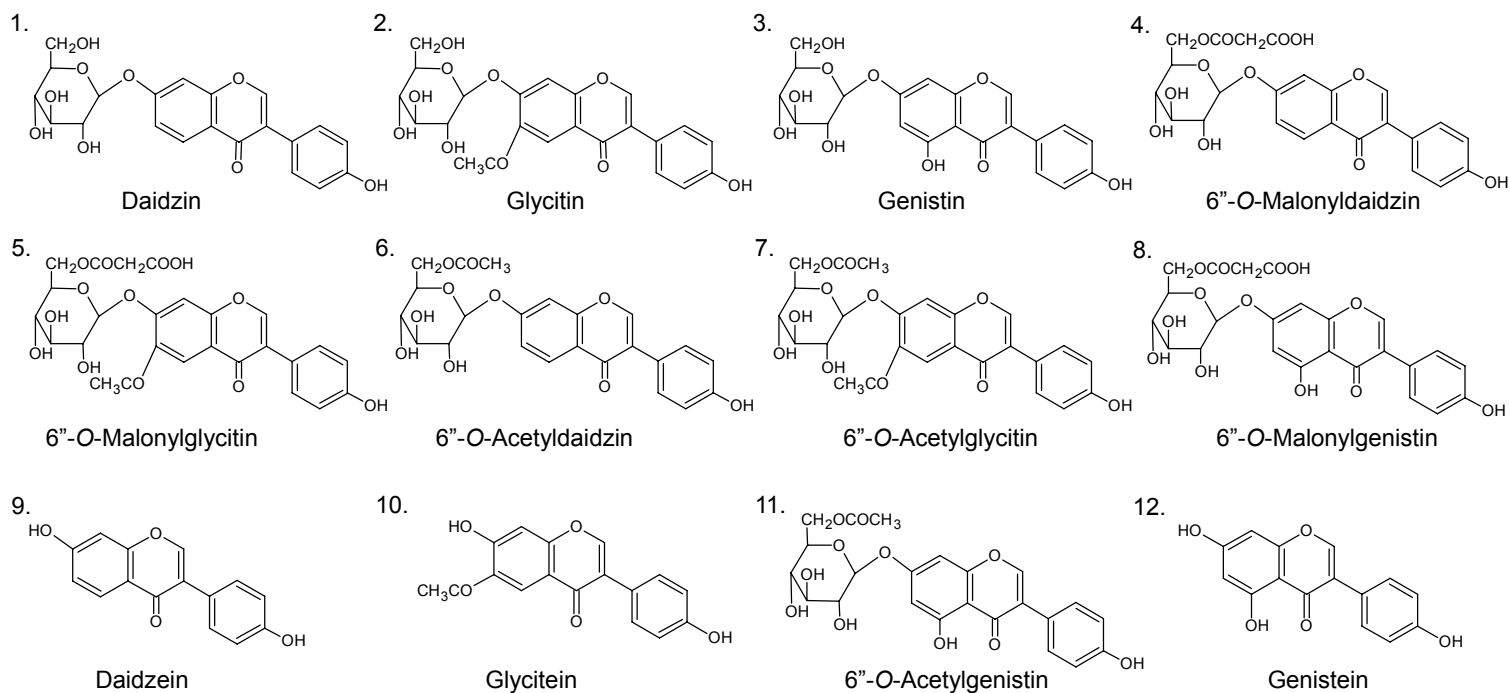


大豆イソフラボン
Soy isoflavones

U100712A

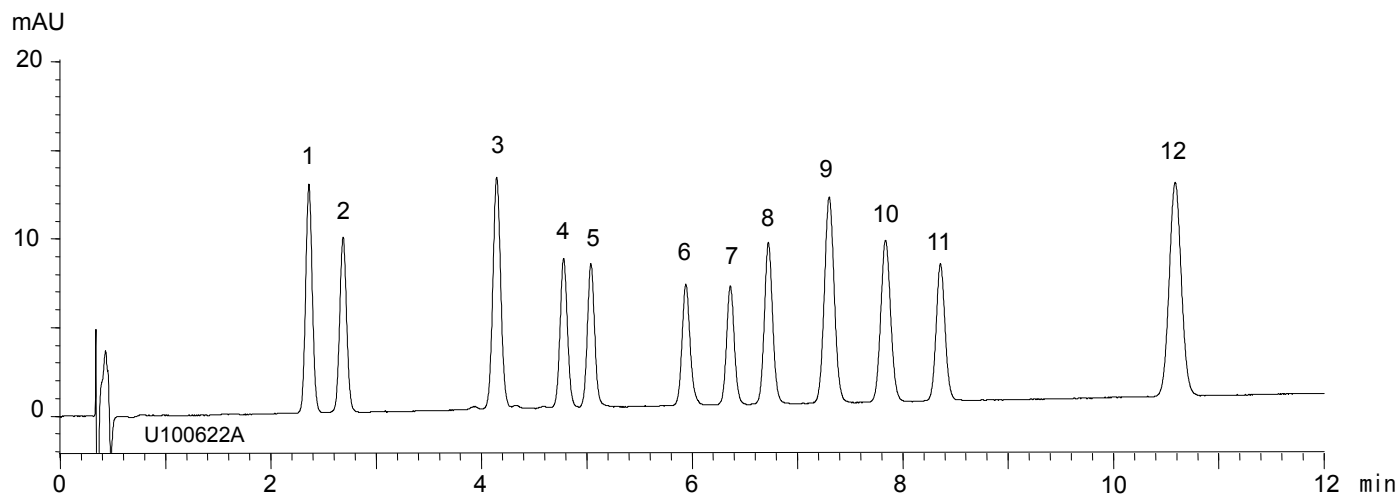

Column : YMC-Triart C18 (3 μ m, 12 nm)
50 X 2.0 mmI.D.

Eluent : A) acetonitrile/water/HCOOH (10/90/0.1)
B) acetonitrile/water/HCOOH (60/40/0.1)
5-40%B (0-12 min)

Flow rate : 0.4 mL/min

Temperature : 25°C

Detection : UV at 254 nm

Injection : 2 μ L

Sample : 0.01 mg/mL