
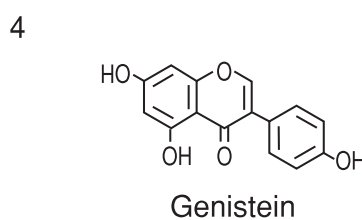
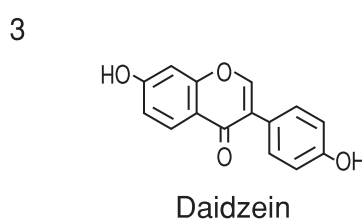
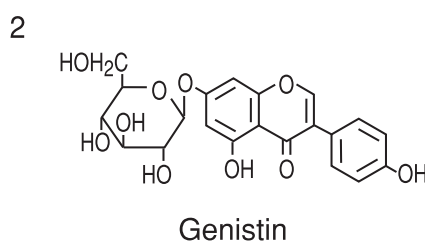
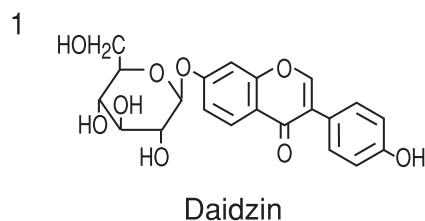


市販みそ中のイソフラボン

Isoflavones from miso, fermented soybean paste

F001025B

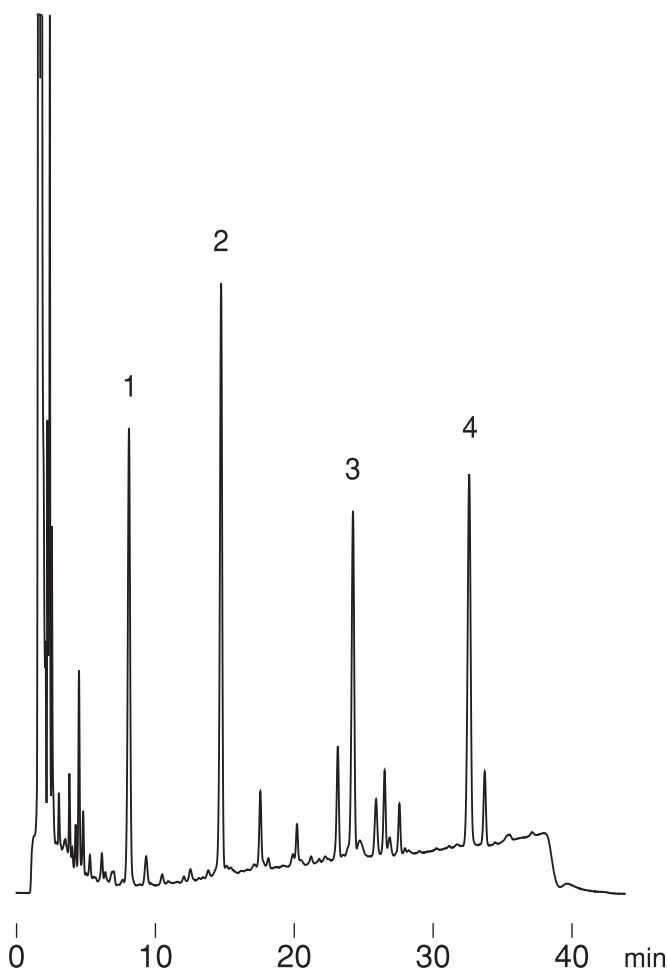

Column : Hydrosphere C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.

Eluent : A)water/acetic acid (100/1)
B)acetonitrile/water/acetic acid (50/50/1)
30%B(0-5min), 30-70%B(5-35min)

Flow rate : 1.0 mL/min

Temperature : 37 $^{\circ}$ C

Detection : UV at 254 nm, 0.08AUFS

Injection : 7 μ L

Sample : methanol extract of miso