



酸化防止剤
Antioxidants

A970521B

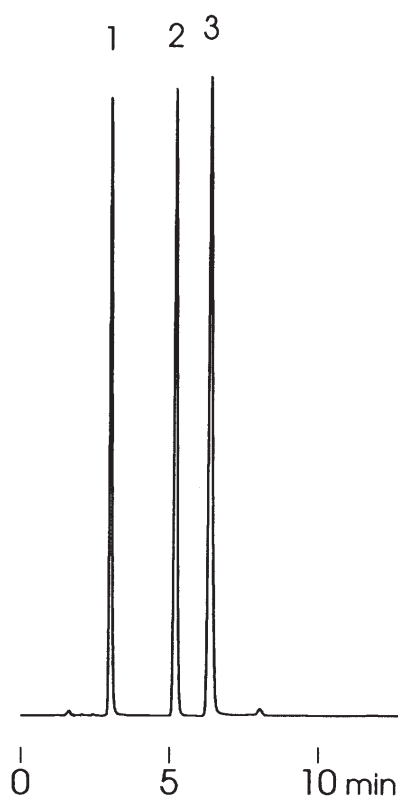


1

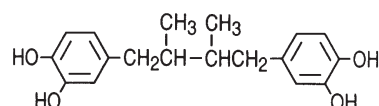


Nordihydroguaiaretic acid(NDGA)

2

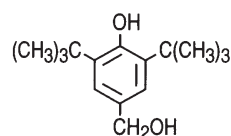

4-Hydroxymethyl-2,6-di-tert-butylphenol (HMBP)

3

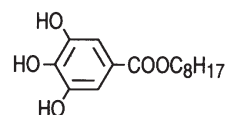

n-Octyl gallate(OG)

Column	: YMC-Pack Pro C18 (5 μm, 12nm) 150 × 4.6 mm I.D.
Eluent	: acetonitrile/methanol/5% acetic acid (35/35/30)
Flow rate	: 1.0 mL/min
Temperature	: 30 °C
Detection	: UV at 280 nm, 0.32 AUFS
Injection	: 10 μL (0.1~0.6 mg/mL)