


多環芳香族炭化水素
Polyaromatic hydrocarbons

K930918A

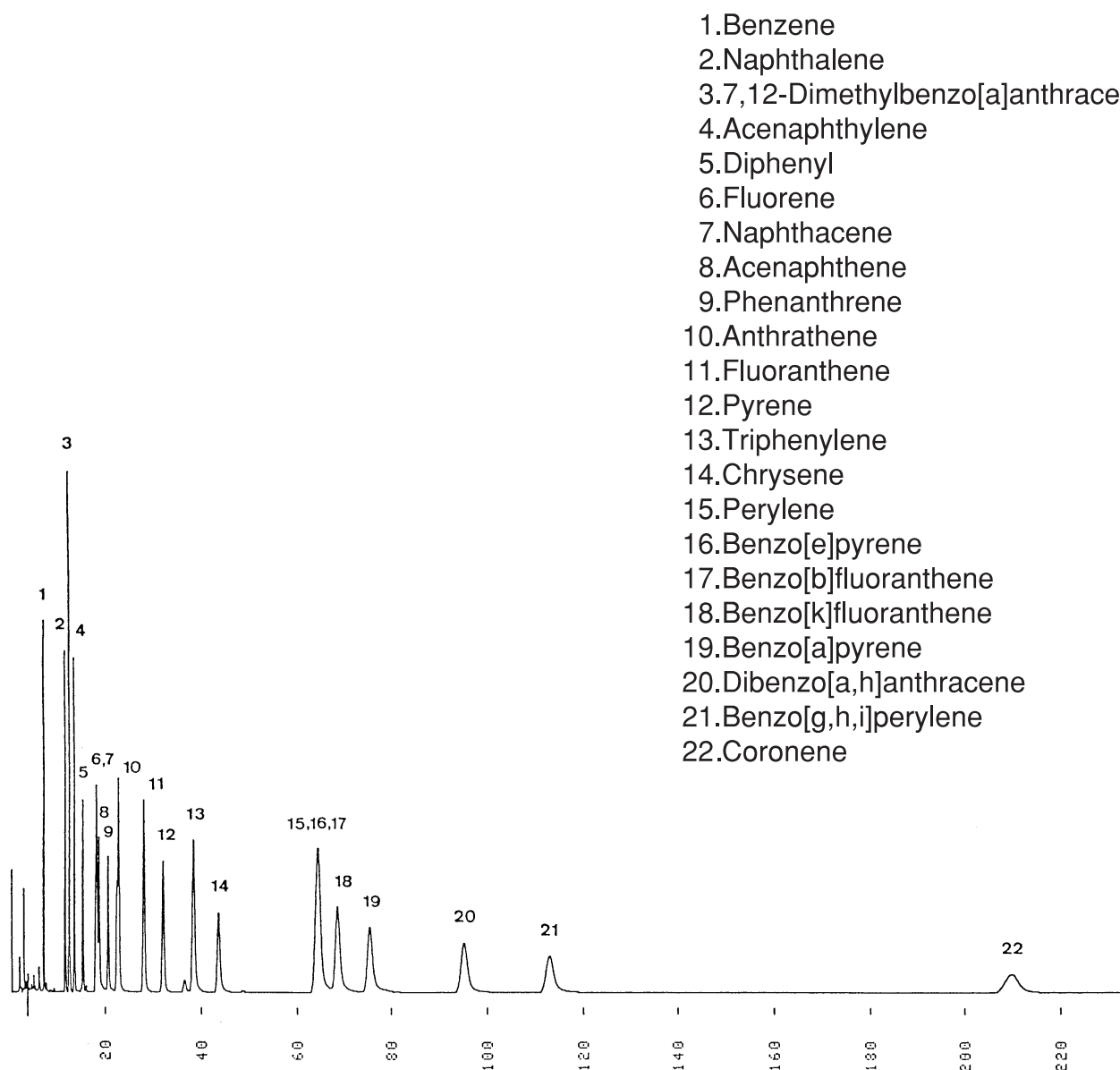

1. Benzene
2. Naphthalene
3. 7,12-Dimethylbenzo[a]anthracene
4. Acenaphthylene
5. Diphenyl
6. Fluorene
7. Naphthacene
8. Acenaphthene
9. Phenanthrene
10. Anthrathene
11. Fluoranthene
12. Pyrene
13. Triphenylene
14. Chrysene
15. Perylene
16. Benzo[e]pyrene
17. Benzo[b]fluoranthene
18. Benzo[k]fluoranthene
19. Benzo[a]pyrene
20. Dibenzo[a,h]anthracene
21. Benzo[g,h,i]perylene
22. Coronene

Column : J'sphere ODS-H80 (4 μ m, 8nm)
250 \times 4.6 mm I.D.
Eluent : acetonitrile/water (65/35)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.08 AUFS
Injection : 7 μ L