


多環芳香族炭化水素
Polycyclic aromatic hydrocarbons

K930918B

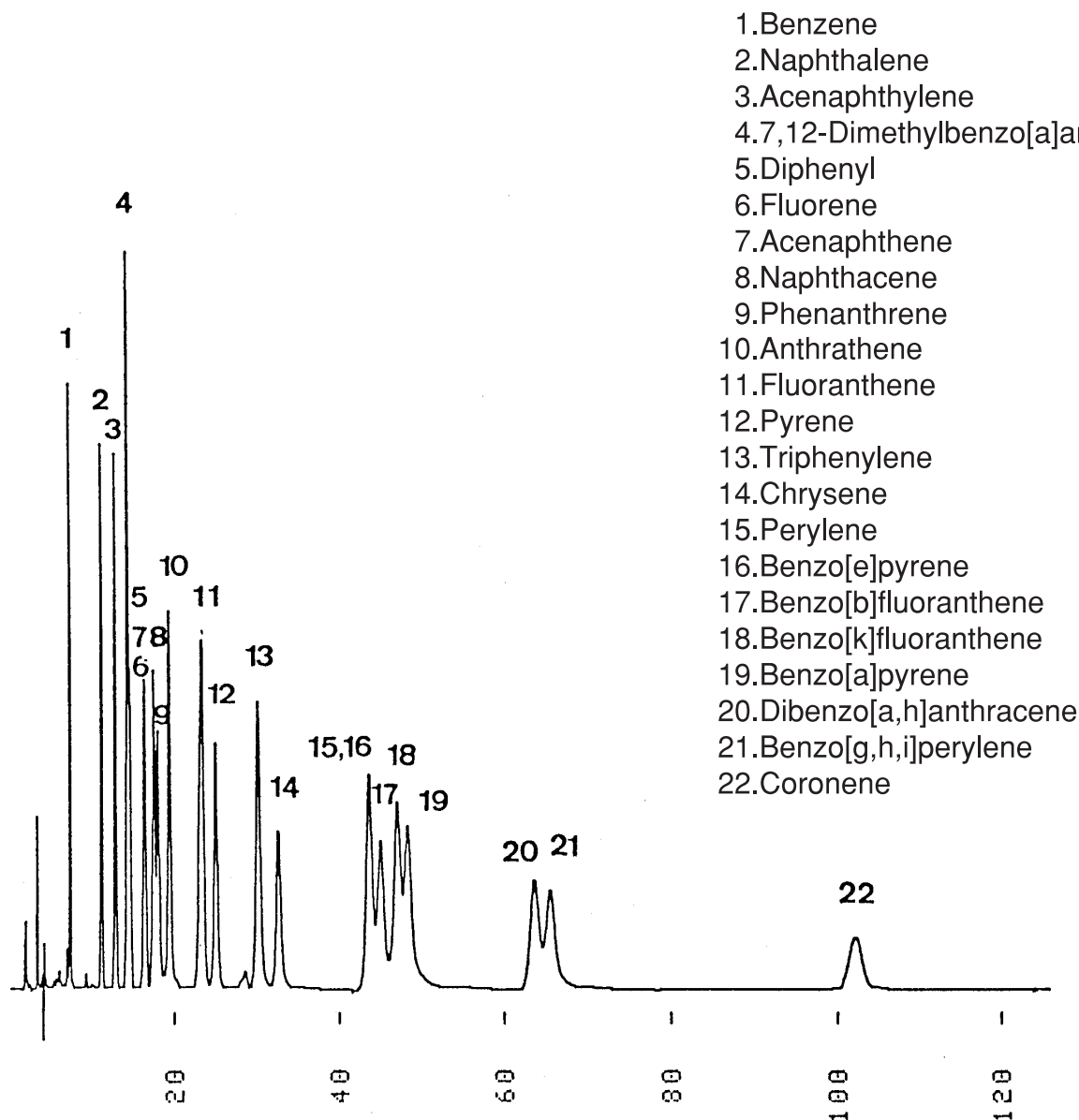

1. Benzene
2. Naphthalene
3. Acenaphthylene
4. 7,12-Dimethylbenzo[a]anthracene
5. Diphenyl
6. Fluorene
7. Acenaphthene
8. Naphthacene
9. Phenanthrene
10. Anthracene
11. Fluoranthene
12. Pyrene
13. Triphenylene
14. Chrysene
15. Perylene
16. Benzo[e]pyrene
17. Benzo[b]fluoranthene
18. Benzo[k]fluoranthene
19. Benzo[a]pyrene
20. Dibenzo[a,h]anthracene
21. Benzo[g,h,i]perylene
22. Coronene

Column : J'sphere ODS-M80 (4 μ m, 8nm)
250 \times 4.6 mm I.D.
Eluent : acetonitrile/water (65/35)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.08 AUFS
Injection : 7 μ L