

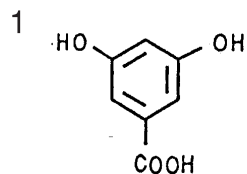


レゾルシン酸異性体  
Resorcylic acid isomers


O940607C

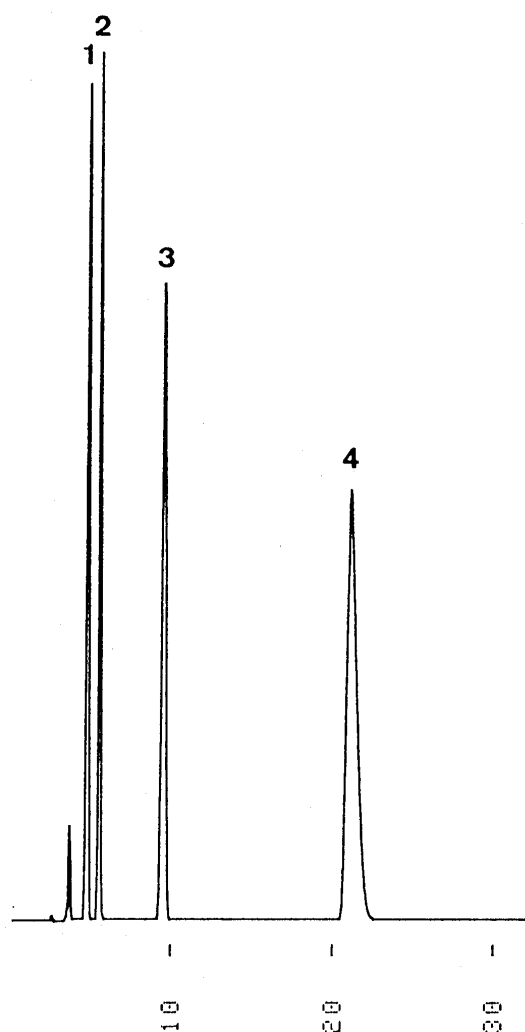


$\alpha$ -Resorcylic acid

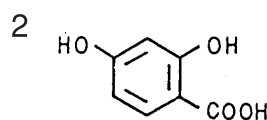

$\beta$ -Resorcylic acid

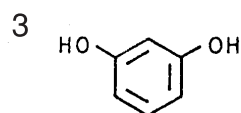

Resorcinol

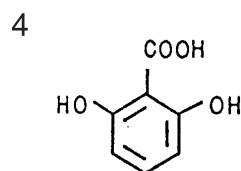

$\gamma$ -Resorcylic acid

Column : YMC-Pack C8 (5  $\mu$ m, 12nm)  
150  $\times$  4.6 mm I.D.  
Eluent : methanol/100mM KH<sub>2</sub>PO<sub>4</sub> (10/90)  
Flow rate : 0.8 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 285 nm, 0.16 AUFS  
Injection : 10  $\mu$ L (0.05~0.4 mg/mL)