
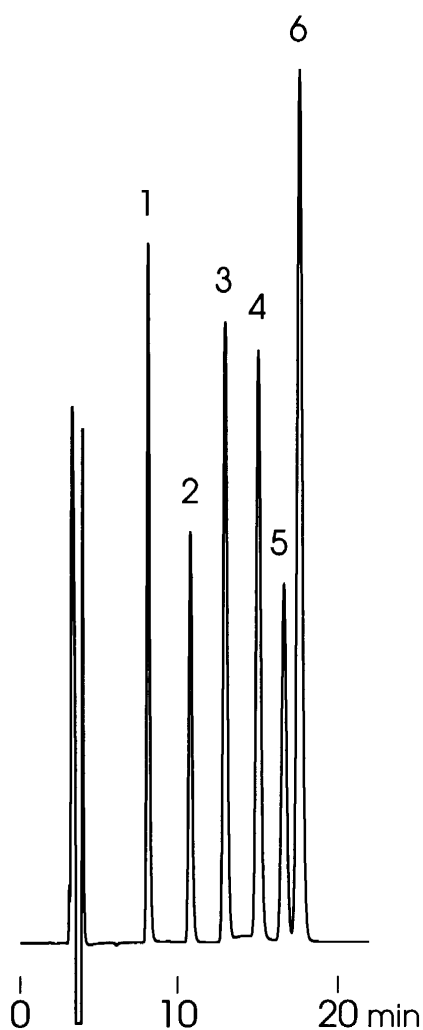
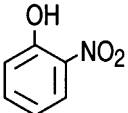
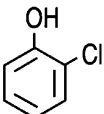
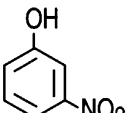
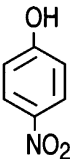



フェノール誘導体
Phenol derivatives

O950614F

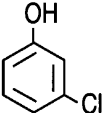

- 1 
o-Nitrophenol
- 2 
o-Chlorophenol
- 3 
m-Nitrophenol
- 4 
p-Nitrophenol
- 5 
m-Chlorophenol
- 6 
p-Chlorophenol

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : tetrahydrofuran/water/acetic acid (35/65/1)
Flow rate : 0.5 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 240 nm, 0.13 AUFS
Injection : 4 μ L (0.08 mg/mL)