


フェノール誘導体
Phenol derivatives

O950614C

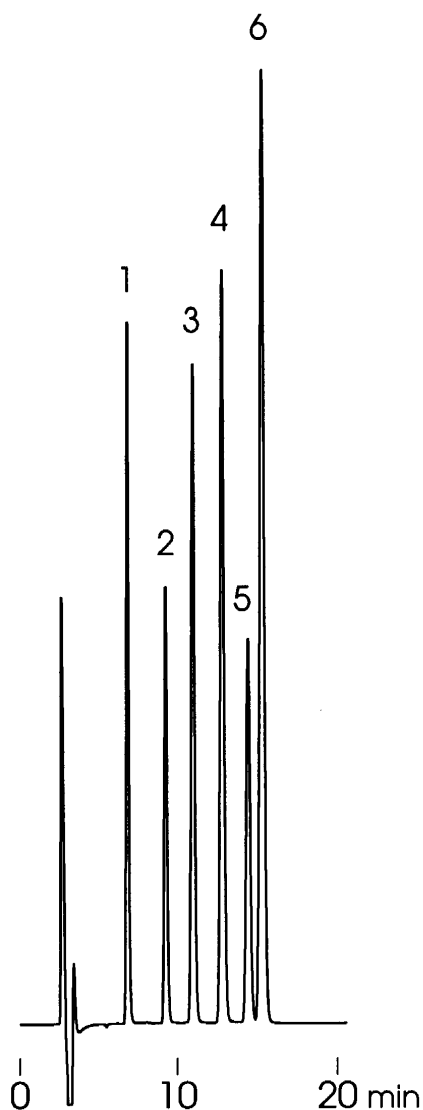

- 1
o-Nitrophenol
- 2
o-Chlorophenol
- 3
m-Nitrophenol
- 4
p-Nitrophenol
- 5
m-Chlorophenol
- 6
p-Chlorophenol

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : tetrahydrofuran/water/acetic acid (35/65/1)
Flow rate : 0.5 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 240 nm, 0.13 AUFS
Injection : 4 μ L (0.08 mg/mL)