



ターフェニル異性体
Terphenyl isomers

J020719E

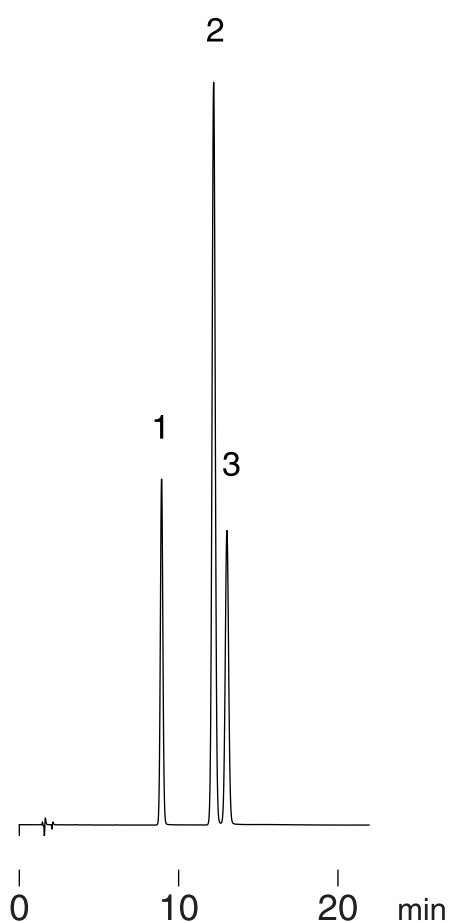


1.

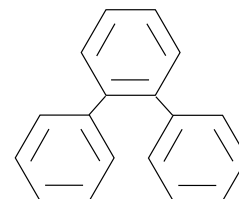


o-Terphenyl

2.

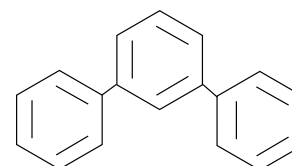

m-Terphenyl

3.

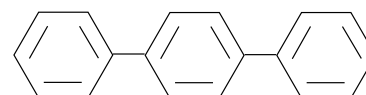

p-Terphenyl

| | |
|-------------|--|
| Column | : YMC-Pack Pro C18 RS (5 μ m, 8nm) 150 \times 4.6 mm I.D. |
| Eluent | : methanol/water (85/15) |
| Flow rate | : 1.0 mL/min |
| Temperature | : 37 $^{\circ}$ C |
| Detection | : UV at 254 nm, 0.16 AUFS |
| Injection | : 10 μ L (0.02 mg/mL) |