

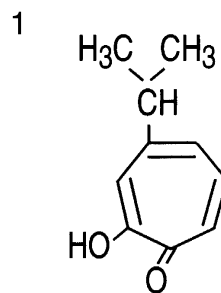


ヒノキチオール
Hinokitiol

O960507F

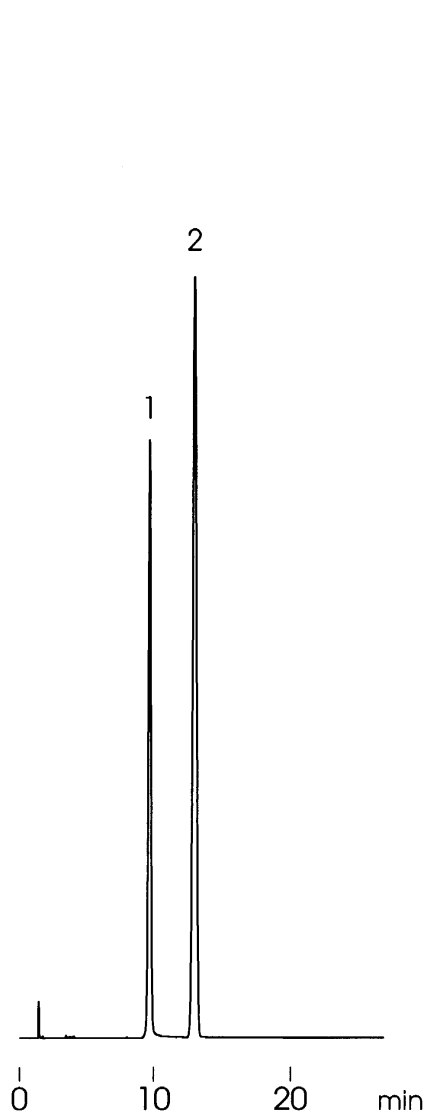

Hinokitiol

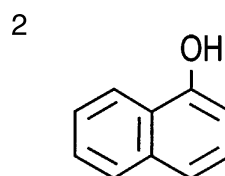

1-Naphthol(I.S.)

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 20mM NH₄H₂PO₄ containing 1mM EDTA \cdot 2Na/acetonitrile (65/35)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 320 nm, 0.16 AUFS
Injection : 10 μ L (0.08, 0.25 mg/mL)