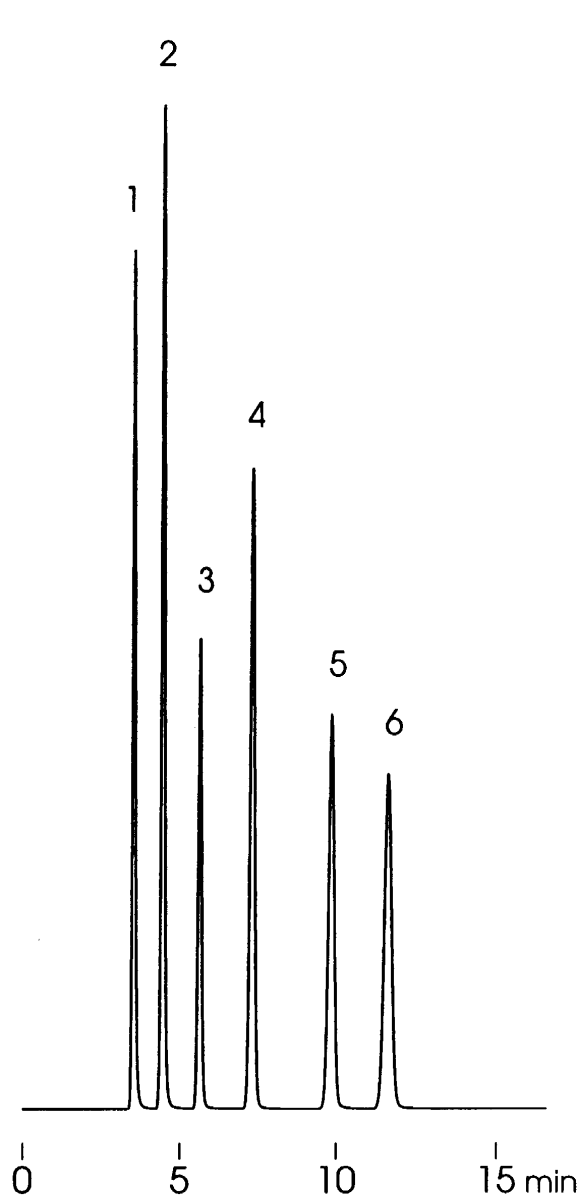
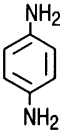

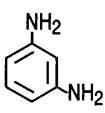
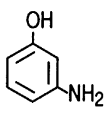
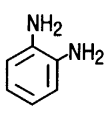
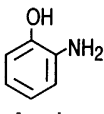


芳香族アミン
Aromatic amines

S961007C



- 1 
p-Phenylenediamine
- 2 
p-Aminophenol
- 3 
m-Phenylenediamine
- 4 
m-Aminophenol
- 5 
o-Phenylenediamine
- 6 
o-Aminophenol

Column : YMC-Pack Pro C18 (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 50mM NH₄H₂PO₄ - (NH₄)₂HPO₄ (pH 6.7)/methanol (95/5)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.26 AUFS
Injection : 5 μ L (0.125~0.5 mg/mL)