
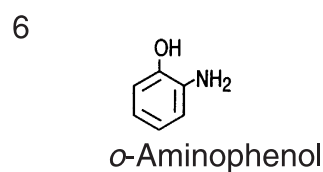
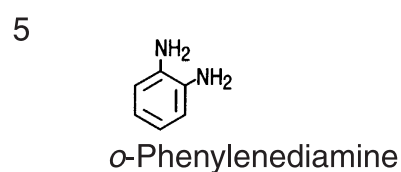
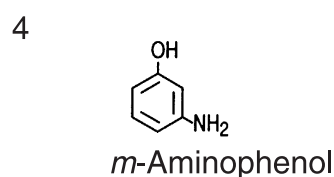
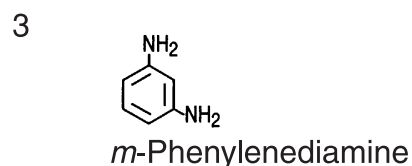
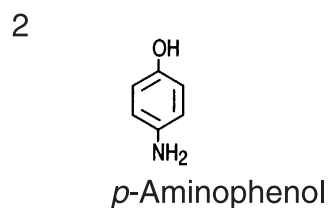
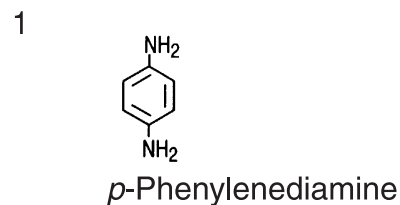


芳香族アミン
Aromatic amines

S961004E

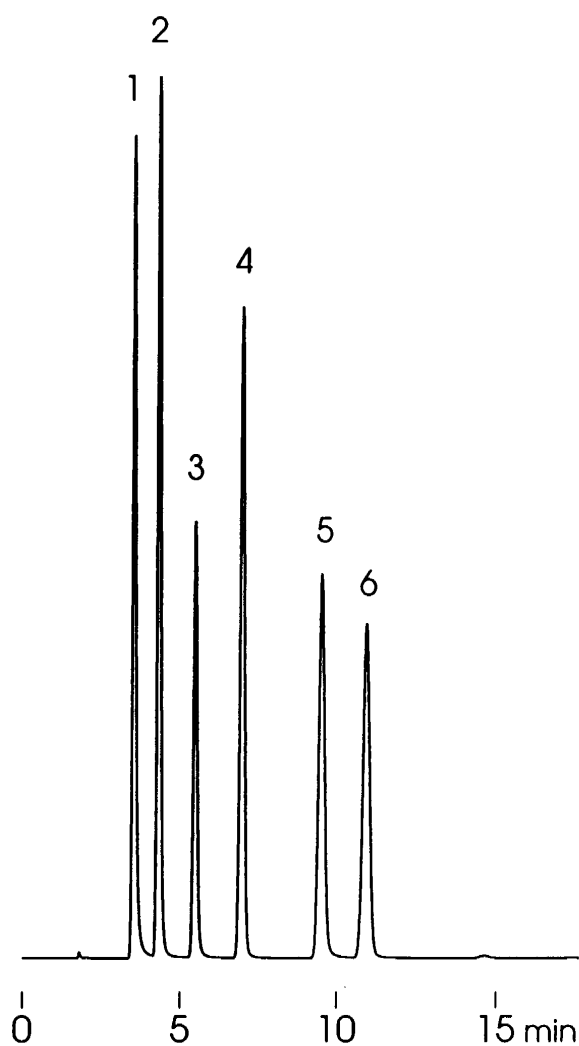

Column : YMC-Pack ODS-AM (5 μ m, 12nm)
150 \times 4.6 mm I.D.
Eluent : 50mM NH₄H₂PO₄ - (NH₄)₂HPO₄ (pH 6.7)/methanol (95/5)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 254 nm, 0.26 AUFS
Injection : 5 μ L (0.125~0.5 mg/mL)