
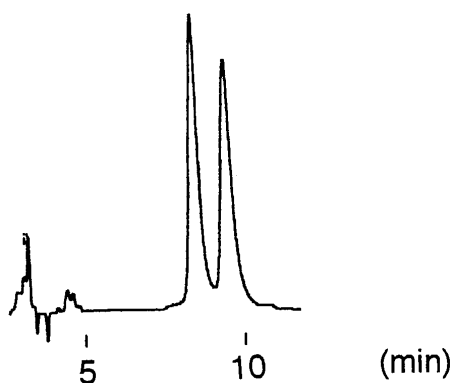


CBZ-フェニルアラニン

CBZ-Phenylalanine (Z-Phe-OH)

N020605N

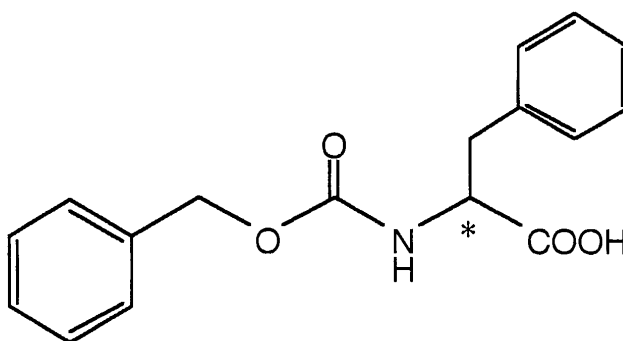

Column	: YMC CHIRAL β -CD BR 250 \times 4.6 mm I.D.
Eluent	: 0.53M AcOH - 0.53M TEA (pH 4.0)/acetonitrile (40/60) (TEA=triethylamine)
Flow rate	: 1.0 mL/min
Temperature	: 25 $^{\circ}$ C
Detection	: UV at 254 nm
Injection	: 5 μ L (1 mg/mL)