
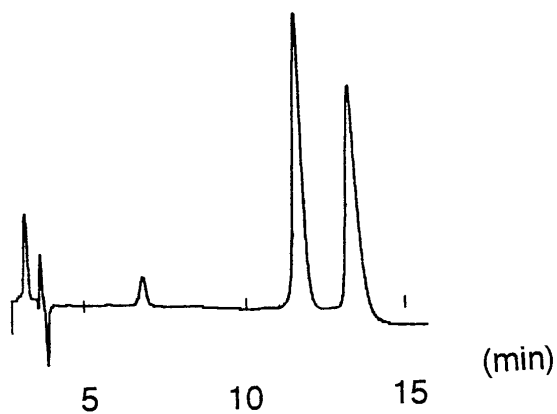


CBZ-フェニルアラニン

CBZ-Phenylalanine (Z-Phe-OH)

N020605T

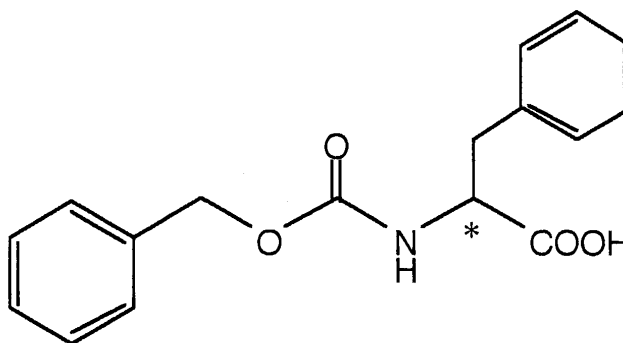

| | |
|-------------|---|
| Column | : YMC CHIRAL γ -CD BR 250 \times 4.6 mm I.D. |
| Eluent | : 0.53M AcOH - 0.53M TEA (pH 4.0)/acetonitrile (60/40) (TEA=triethylamine) |
| Flow rate | : 1.0 mL/min |
| Temperature | : 25 $^{\circ}$ C |
| Detection | : UV at 254 nm |
| Injection | : 5 μ L (1 mg/mL) |