

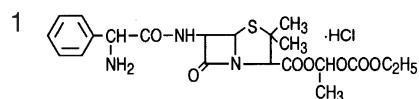


ペニシリン系抗生物質
Penicillin antibiotics


O950908A

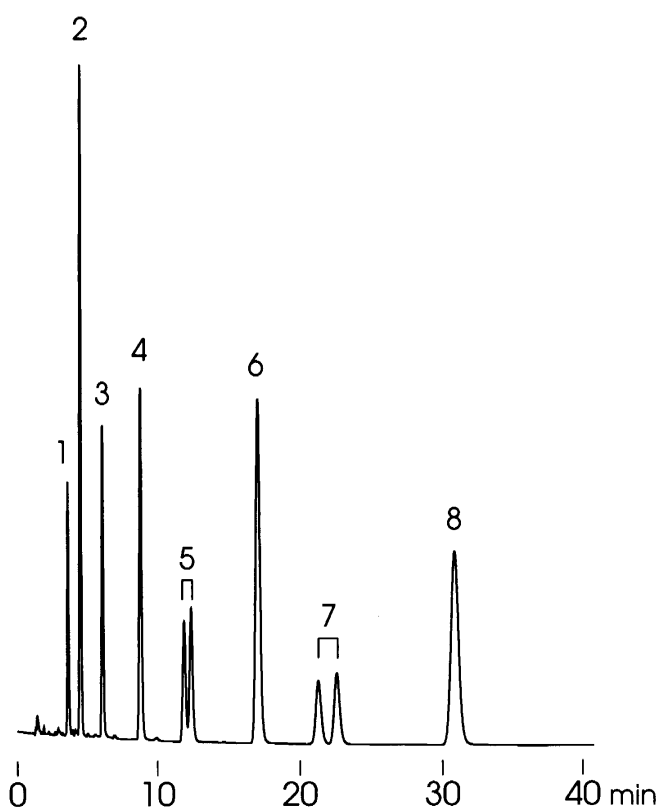


Bacampicillin hydrochloride

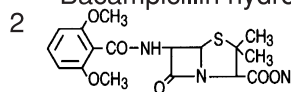


Methicillin sodium

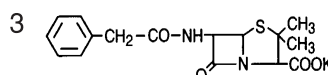


Penicillin G potassium
(Benzylpenicillin potassium)

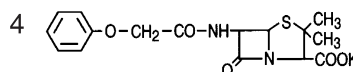


Penicillin V potassium
(Phenoxymethylpenicillin potassium)

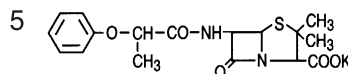


Phenethicillin potassium
(Phenoxyethylpenicillin potassium)

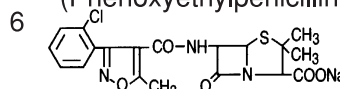

Cloxacillin sodium

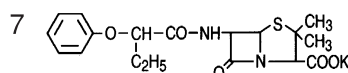

Propicillin potassium

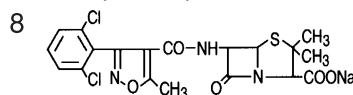

Dicloxacillin sodium

Column : J'sphere ODS-H80 (4 μ m, 8nm)
150 \times 4.6 mm I.D.
Eluent : 50mM NaH₂PO₄ - H₃PO₄ (pH 2.3)/acetonitrile (67/33)
Flow rate : 1.0 mL/min
Temperature : 37 $^{\circ}$ C
Detection : UV at 220 nm, 0.26 AUFS
Injection : 8 μ L (0.13 mg/mL)