


クリンダマイシン塩酸塩（米国薬局方記載条件）

Clindamycin hydrochloride (The United States Pharmacopoeia)


F100416A

(A) Standard solution

(0.05 mg/mL Lincomycin HCl, 0.1 mg/mL Clindamycin HCl)

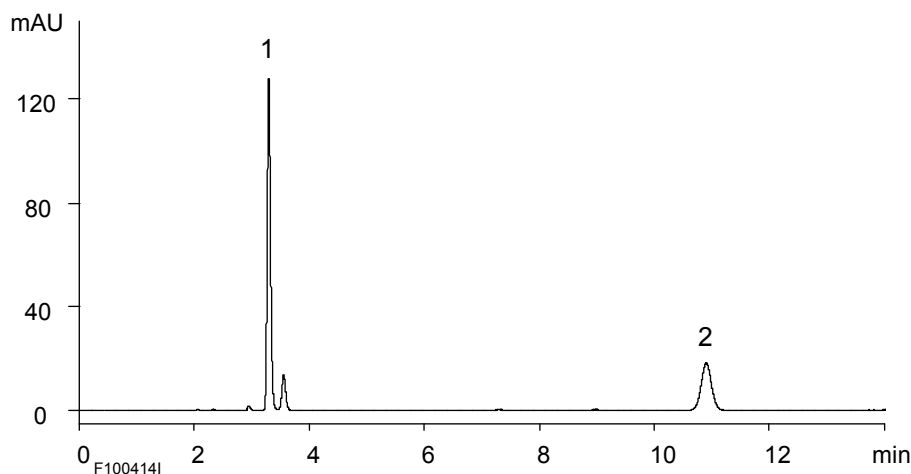


1.

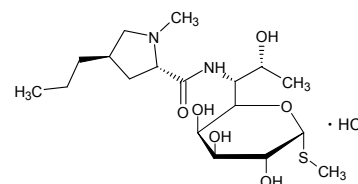


Lincomycin hydrochloride

(B) Test solution

(5.0 mg/mL Clindamycin HCl)

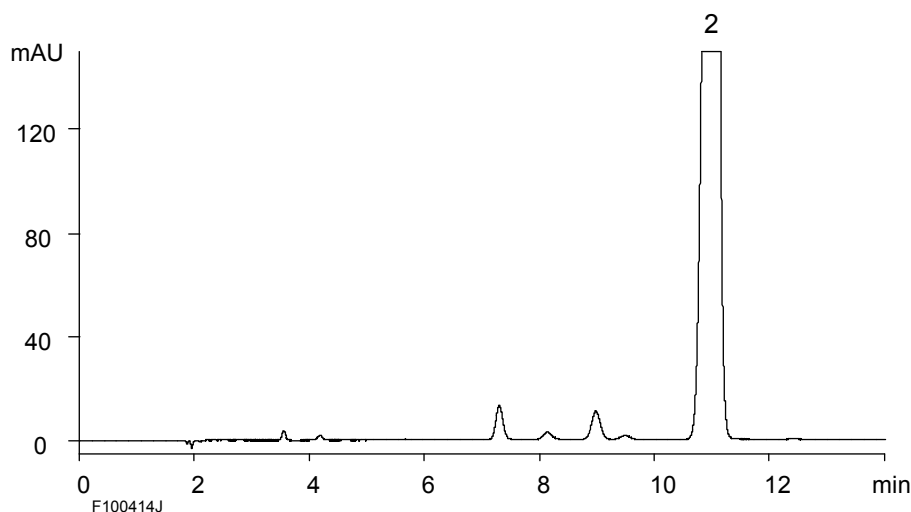

2.

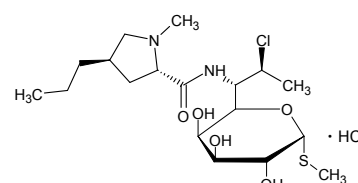

Clindamycin hydrochloride

Column	: YMC-Triart C18 (5 μ m, 12 nm) 250 X 4.6 mm I.D.
Eluent	: 50 mM KH_2PO_4 (pH 7.5 adjusted by 8 M KOH)/acetonitrile (55/45)
Flow rate	: 1.0 mL/min
Temperature	: 25°C
Detection	: UV at 210 nm
Injection	: (A) 10 μ L (B) 10 μ L