




ビピリジニウム系除草剤(ジクワット、パラコート)
Bipyridinium herbicides (diquat, paraquat)

R091201B

1.

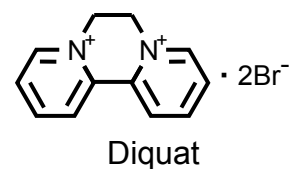

2.

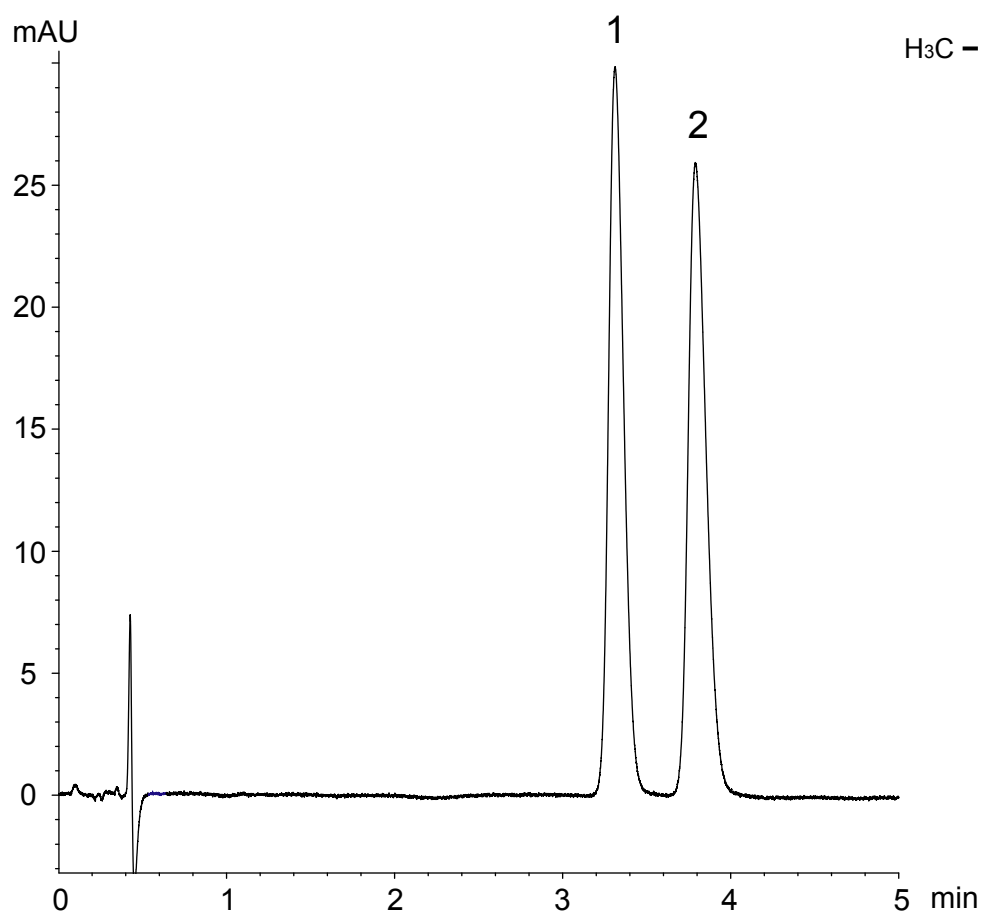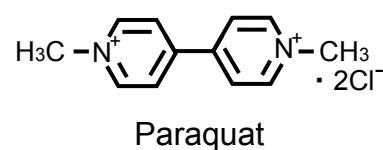

Column : YMC-Triart C18 (3 μ m, 12 nm)
50 X 2.0 mmI.D.
Eluent : water/acetonitrile/HFBA* (95/5/0.1)
Flow rate : 0.4 mL/min
Temperature : 37°C
Detection : UV at 290 nm
Injection : 1 μ L (0.1 mg/mL)

*heptafluorobutyric acid