


除草剂  
Herbicides

O941018B

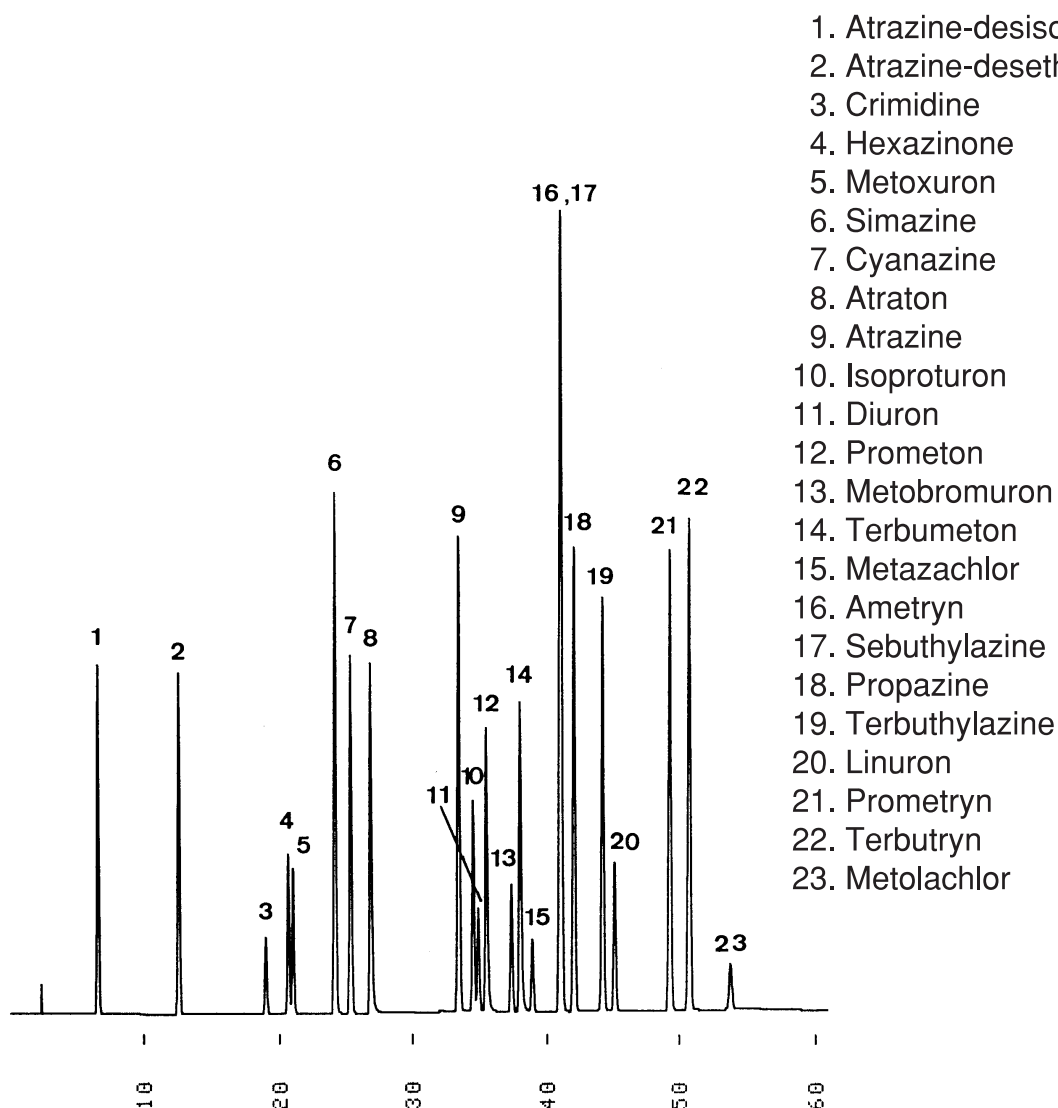

1. Atrazine-desisopropyl
2. Atrazine-desethyl
3. Crimidine
4. Hexazinone
5. Metoxuron
6. Simazine
7. Cyanazine
8. Atraton
9. Atrazine
10. Isoproturon
11. Diuron
12. Prometon
13. Metobromuron
14. Terbumeton
15. Metazachlor
16. Ametryn
17. Sebuthylazine
18. Propazine
19. Terbuthylazine
20. Linuron
21. Prometryn
22. Terbutryn
23. Metolachlor

Column : J'sphere ODS-H80 (4  $\mu$ m, 8nm)  
250  $\times$  4.6 mm I.D.

Eluent : A)acetonitrile/50mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (20/80)  
B)acetonitrile/50mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (65/35)  
0-100%B (0-60min, linear)

Flow rate : 1.0 mL/min

Temperature : 37 °C

Detection : UV at 230 nm, 0.16 AUFS

Injection : 5  $\mu$ L (0.04 mg/mL)