
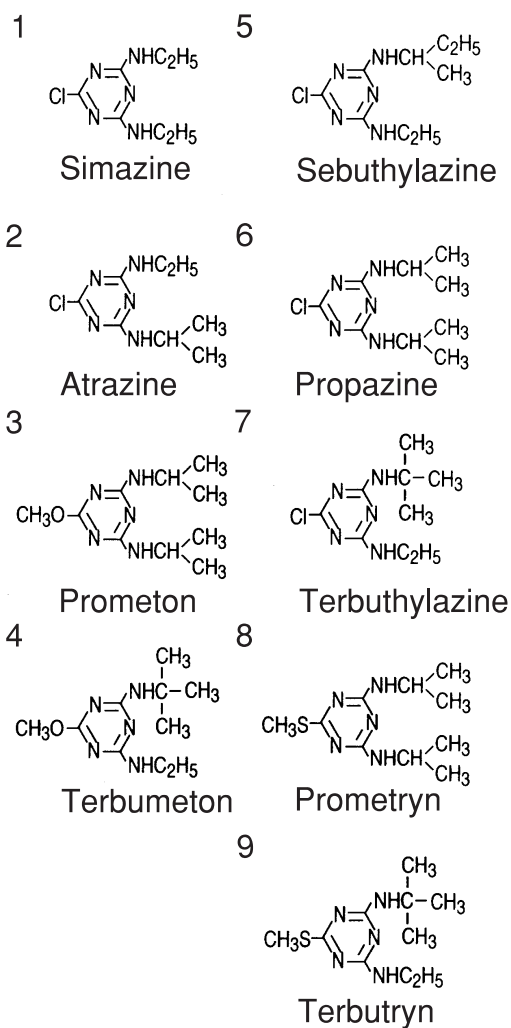


トリアジン系除草剤  
Triazine herbicides

A970210B

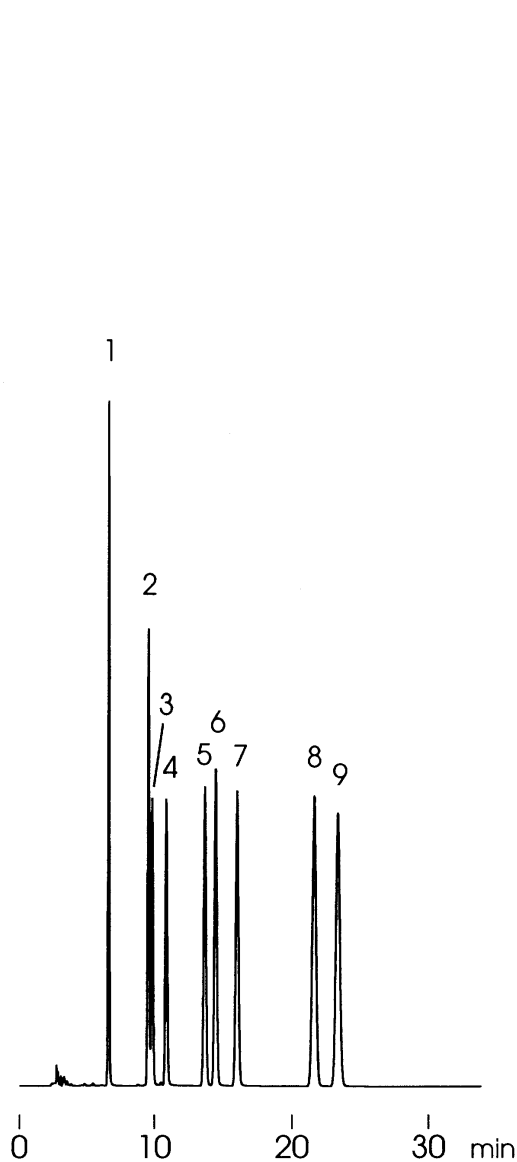

Column : YMC-Pack Pro C18 (5  $\mu$ m, 12nm)  
250  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/50mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (45/55)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.32 AUFS  
Injection : 5  $\mu$ L (0.05 mg/mL)