


トリアジン系除草剤  
Triazine herbicides

O941102F

1. Cyanazine
2. Simazine
3. Atraton
4. Atrazine
5. Prometon
6. Terbumeton
7. Sebuthylazine
8. Propazine
9. Terbutylazine
10. Prometryn
11. Terbutryn

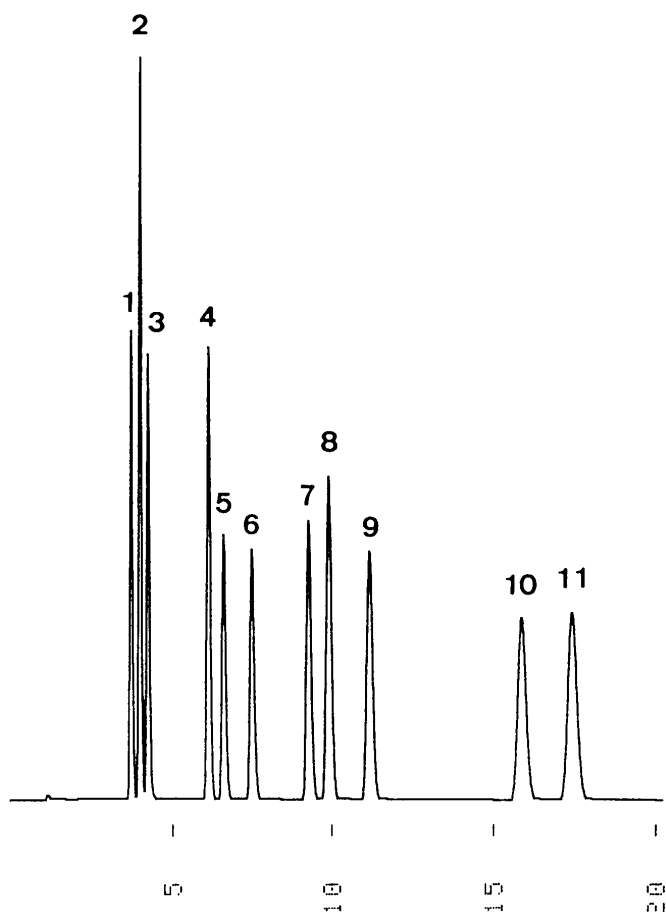

Column : J'sphere ODS-H80 (4  $\mu$  m, 8nm)  
150  $\times$  4.6 mm I.D.  
Eluent : acetonitrile/50mM NH<sub>4</sub>H<sub>2</sub>PO<sub>4</sub> (45/55)  
Flow rate : 1.0 mL/min  
Temperature : 37  $^{\circ}$ C  
Detection : UV at 230 nm, 0.32 AUFS  
Injection : 5  $\mu$ L (0.05 mg/mL)