

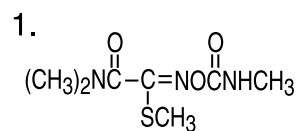


カルバメート系殺虫剤
Carbamate insecticides


J010220D

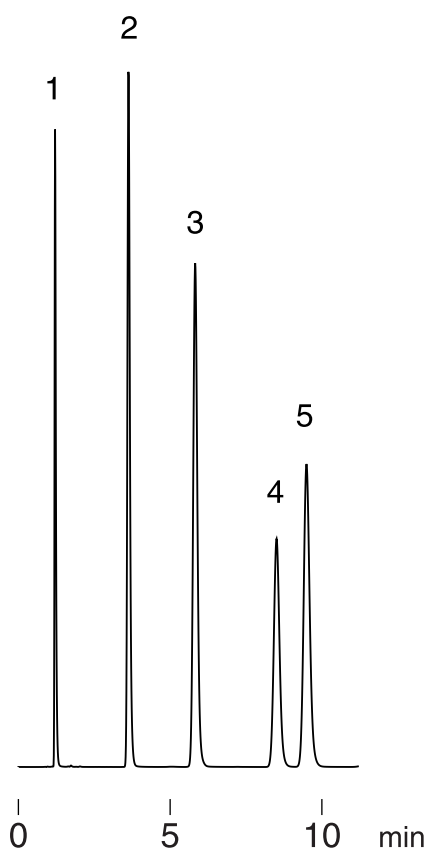


Oxamyl

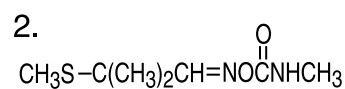


Aldicarb

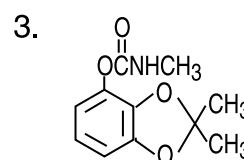

Bendiocarb

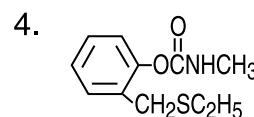

Ethiofencarb

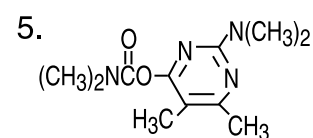

Pirimicarb

Column : Hydrosphere C18 (3 μ m, 12nm)
75 \times 4.6 mm I.D.
Eluent : methanol / 20mM KH₂PO₄ (45/55)
Flow rate : 1.0 mL/min
Temperature : 30 $^{\circ}$ C
Detection : UV at 240 nm, 0.64 AUFS
Injection : 5 μ L (0.1~1.5 mg/mL)