

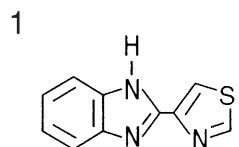


防かび剤  
Fungicides


A970217A

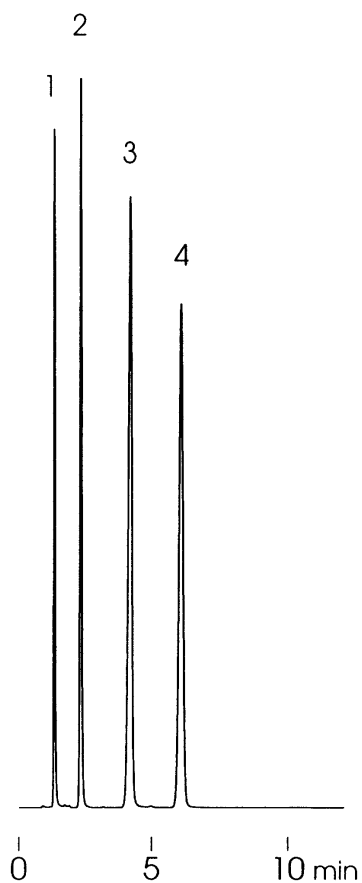


Thiabendazole

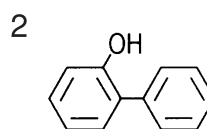

o-Phenylphenol

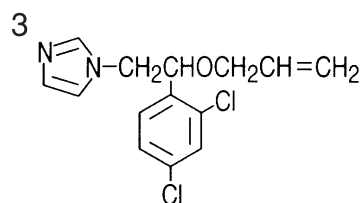

Imazalil

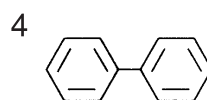

Diphenyl

Column	: YMC-Pack Pro C18 (5 $\mu$ m, 12nm) 75 $\times$ 4.6 mm I.D.
Eluent	: methanol/20mM NH <sub>4</sub> H <sub>2</sub> PO <sub>4</sub> (70/30)
Flow rate	: 1.0 mL/min
Temperature	: 37 $^{\circ}$ C
Detection	: UV at 230 nm, 0.32 AUFS
Injection	: 10 $\mu$ L (0.02~0.20 mg/mL)